

Choose Excellence,
Expect Value

SharpPoint[™] **PLUS**

Surgical Sutures

Suture Catalog



Choose Excellence, Expect Value

Surgical Specialties Corporation is the world's largest independent suture manufacturer. While headquartered in Westwood, Massachusetts, SSC is a global company with development and manufacturing facilities in the United States, Germany, Mexico, China, and the United Kingdom.

We design and manufacture high performance wound closure products including Quill®, Sharpoint™, and LOOK™ brand sutures.

We are pleased to offer the new Sharpoint™ PLUS suture portfolio featuring German engineered needles.

Our philosophy is that performance doesn't have to be expensive.



Table of Contents

| | |
|--|-------|
| Suture Material Selection | 3-4 |
| Advanced Needle Technology | 5 |
| Needle Anatomy and Terminology | 6 |
| Needle Points and Cutting Edges | 7 |
| Needle Charts..... | 8 |
| Absorbable Sutures | |
| Express Gut..... | 9 |
| Plain Gut..... | 9-10 |
| Chromic Gut | 11-12 |
| PolySyn FA™ (polyglycolic acid) | 13 |
| Monoderm™ (glycolide and ε-caprolactone) | 14-15 |
| PolySyn™ (polyglycolic acid)..... | 16-19 |
| PDO™ (polydioxanone) | 20-21 |
| Non-Absorbable Sutures | |
| Silk..... | 22-23 |
| Nylon™ | 24-25 |
| Polypropylene..... | 26-29 |
| Polyviolen™ XL (polyester)-Coated..... | 30-31 |
| Cross Reference Guide | |
| Ethicon Cross Reference..... | 32-41 |
| Covidien Cross Reference..... | 42-46 |



Guide to Suture Materials

Absorbable

| Material | Material Description | BSR Profile | Mass Absorption |
|--|---|--|---|
| EXPRESS GUT™  Comparable to: Fast Absorbing Gut Suture | Twisted <ul style="list-style-type: none"> Natural Packaged in Wetting Solution Common uses: Skin Closure | In vivo strength retention: 40% at 5 days 0% at 10 days In-Vivo Study performed by NAMSA Data on file at Surgical Specialties Corporation | Essentially complete by 30 days In-Vivo Study performed by NAMSA Data on file at Surgical Specialties Corporation |
| GUT  Comparable to: Surgical Gut Suture - Plain | Twisted <ul style="list-style-type: none"> Plain & Chromic Packaged in Wetting Solution Common uses: General Closure, OB/Gyn, Urology, Bowel Anastomosis | In vivo strength retention: Varies | Varies depending upon type of suture (plain vs. chromic) and site of placement |
| CHROMIC GUT  Comparable to: Surgical Gut Suture - Chromic | Twisted <ul style="list-style-type: none"> Plain & Chromic Packaged in Wetting Solution Common uses: General Closure, OB/Gyn, Urology, Bowel Anastomosis | In vivo strength retention: Varies | Varies depending upon type of suture (plain vs. chromic) and site of placement |
| POLYSYN FA™ (polyglycolic acid)  · VICRYL RAPIDE™ * Velosorb™ Fast | Braided <ul style="list-style-type: none"> Undyed (clear) Coated Common uses: General Closure, OB/Gyn, Urology | In vivo strength retention: 48% at 7 days 0% at 14 days | Essentially complete between 49 and 63 days |
| MONODERM™ (polyglycolic acid and e-caprolactone)  · MONOCRYL™ * Caprosyn™/Biosyn™ | Monofilament <ul style="list-style-type: none"> Undyed (clear) Dyed (violet) Common uses: Skin Closure, OB/Gyn, Urology, Gastrointestinal | In vivo strength retention: 42% to 76% at 7 days 36% to 52% at 14 days | Essentially complete by 90 days |
| POLYSYN™ (polyglycolic acid)  · VICRYL™ * Polysorb™/ Dexon™ | Braided <ul style="list-style-type: none"> Undyed (clear) Dyed (violet) Monofilament (9-0 & 10-0 only) <ul style="list-style-type: none"> Undyed (clear) Dyed (violet) Coated & Uncoated Common uses: Skin Closure, OB/Gyn, Urology, Gastrointestinal | In vivo strength retention: 50% at 14 days 20% at 21 days | Essentially complete between 50 and 90 days |
| PDO™ (polydioxanone)  · PDS II™ * Maxon™ | Monofilament <ul style="list-style-type: none"> Dyed (violet) Common uses: Plastic Surgery, Orthopaedic, OB/Gyn, Gastrointestinal | In vivo strength retention: 80% to 90% at 14 days 60% to 82% at 28 days 47% to 79% at 42 days | Essentially complete between 180 and 220 days |

· Comparable Ethicon™ Material * Comparable Covidien™ Material


Non-Absorbable

| Material | Material Description | Qualities |
|---|---|---|
| SILK  · PERMA-HAND™ Silk * SofSilk™ | Braided <ul style="list-style-type: none"> • Undyed (white) • Dyed (black) • Coated Common uses: General & Skin Closure, Gastrointestinal, Cardiovascular Surgery, Plastic Surgery, Ophthalmology | <ul style="list-style-type: none"> • Excellent handling and tying characteristics |
| NYLON™  · ETHILON™ * Monosof™/ Dermalon™ | Monofilament <ul style="list-style-type: none"> • Undyed (clear) • Dyed (blue and black) • Packaged dry or in wetting solution Common uses: General & Skin Closure, Plastic Surgery, Ophthalmology, Microsurgery, Cardiovascular Surgery, Gastrointestinal | <ul style="list-style-type: none"> • Passes easily through tissue |
| POLYPROPYLENE  · PROLENE™ * Surgipro™/Surgipro™ II | Monofilament <ul style="list-style-type: none"> • Dyed (blue) Common uses: General & Skin Closure, Neurosurgery, Plastic Surgery, Ophthalmology, Microsurgery, Cardiovascular Surgery | <ul style="list-style-type: none"> • Inert • Easy to handle • Passes easily through tissue |
| POLYVIOLENE™ (UNCOATED) <small>(poly (ethylene-terephthate))</small>  · MERSILENE™ * Surgidac™ | Braided <ul style="list-style-type: none"> • Undyed (white) • Dyed (green) • Uncoated Common uses: General & Skin Closure, Plastic Surgery, Ophthalmology, Cardiovascular Surgery, Neurological Procedures | <ul style="list-style-type: none"> • Smooth • Easy to handle • Good knot security |
| POLYVIOLENE™ XL (COATED)  · ETHIBOND™ * Tricon™ | Braided <ul style="list-style-type: none"> • Undyed (white) • Dyed (green) • Coated Common uses: General & Skin Closure, Plastic Surgery, Ophthalmology, Cardiovascular Surgery, Neurological Procedures | <ul style="list-style-type: none"> • Smooth • Easy to handle • Good knot security |

· Comparable Ethicon™ Material * Comparable Covidien™ Material


Advanced Needle Technology

Sharpoint Plus™ needle performance is an output of its design. Materials, geometry, and coatings work together to produce a high performance surgical instrument.

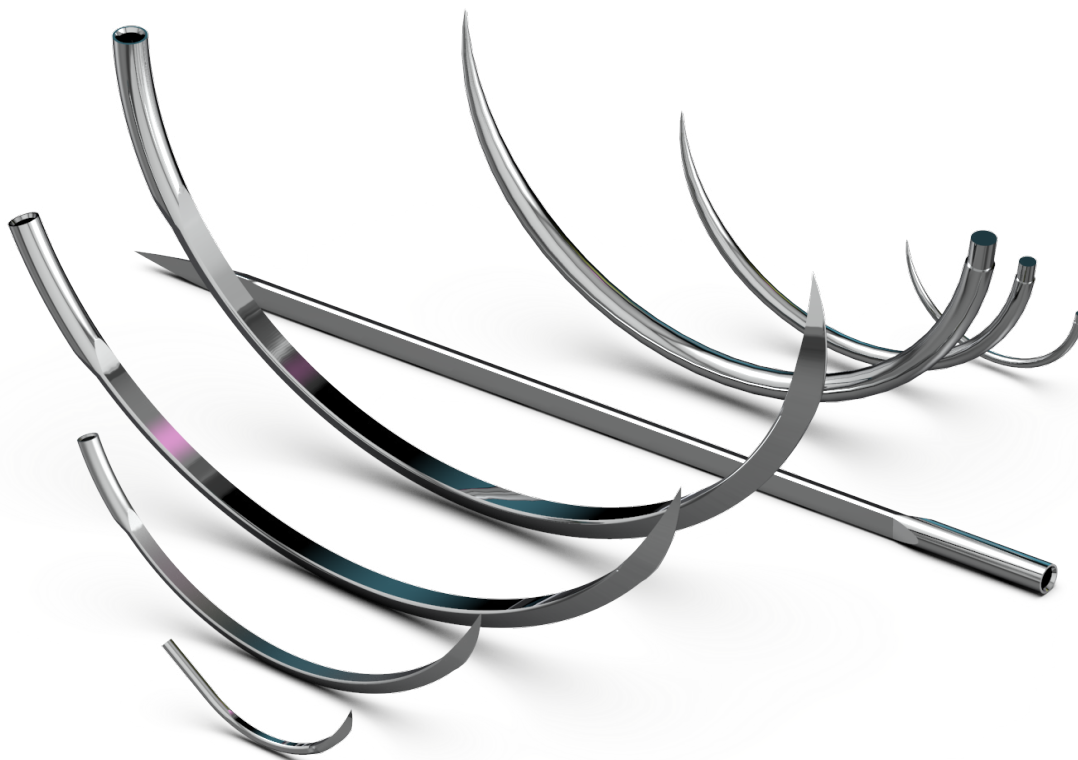


Proprietary Stainless Steel Alloy


Sharpoint Plus™ needles are forged from an alloy of Iron, Nickel, and Chromium. Each element serves a role in delivering strength, ductility, and hardness.

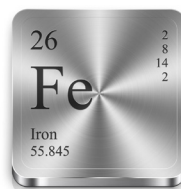

Strength

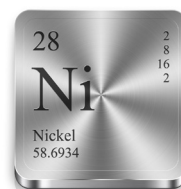

Ductility

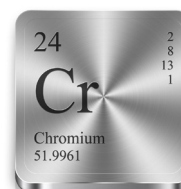

Hardness

Needle Anatomy And Terminology

Needle Terminology

Length: Distance measured along the circumference of the needle, from the point to the swage. It is depicted on the packaging.

Shape: Different surgeries or surgical techniques demand different needle shapes.

Point: The sharpened part of the needle that first penetrates the tissue. The point consists of the area from the sharpened tip of the needle to the widest cross section of the needle body. The point determines ease of penetration and the initial size and shape of the hole.


Swage: Connection point of the suture and the needle. Micro teeth grip the suture material and prevent it from slipping out.

Wire Diameter: Gauge or thickness of the wire originally used to form the needle body.


Body: The part between the point and the swage of the needle.

Included Angle (Curvature): The fraction of the circle included in the curve of the needle. This is depicted on the packaging.


Needle Anatomy

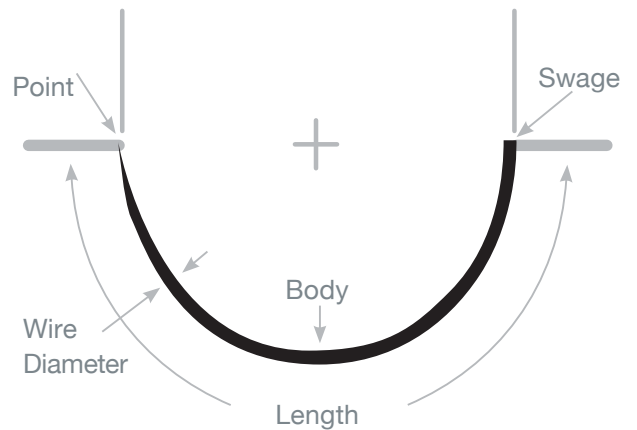


Needle Shapes and General Uses

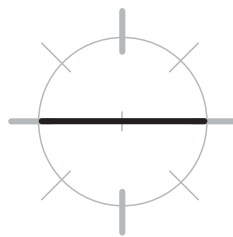


Straight (aka: Keith): Skin Closure

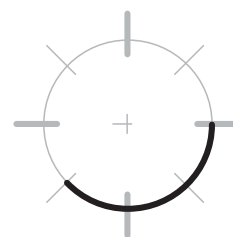

3/8 Circle: Skin Closure, Ophthalmic

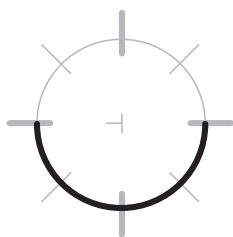



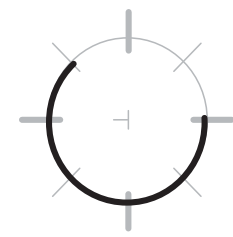


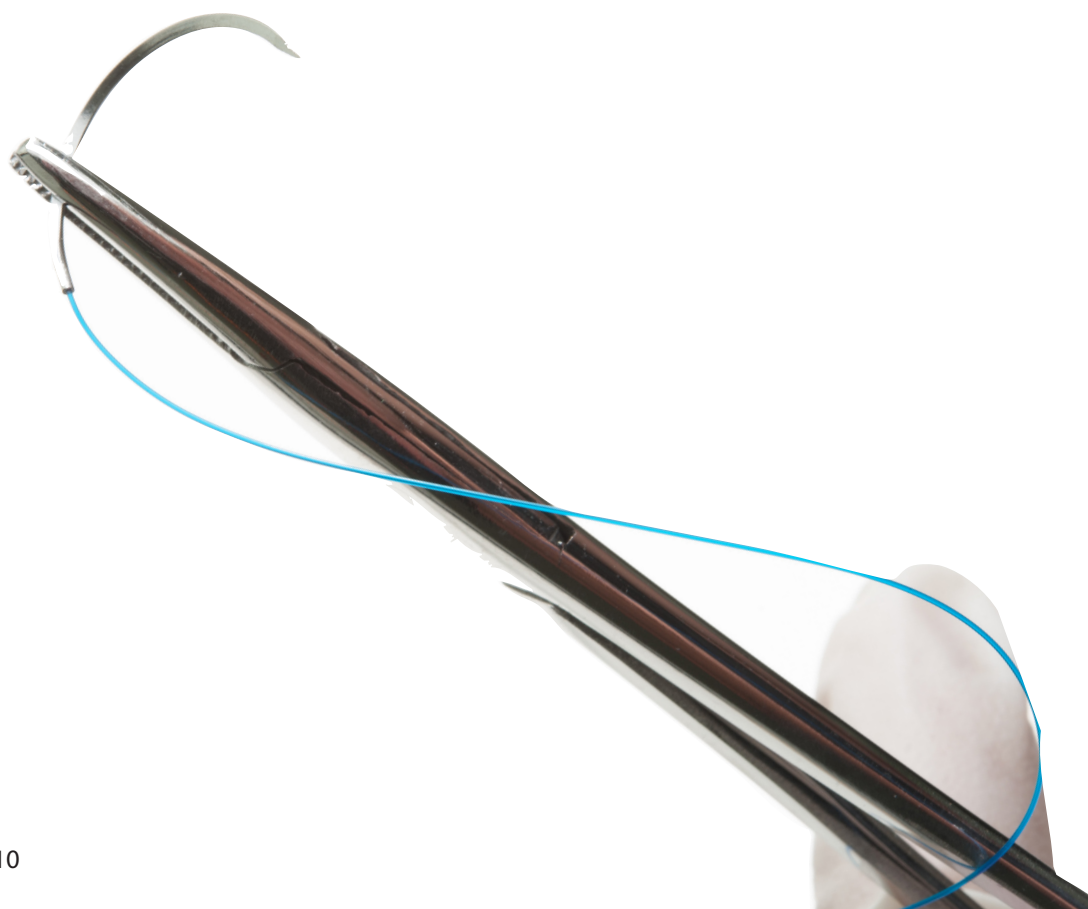
1/2 Circle: Skin Closure, GI, Cardiovascular, Ophthalmic






5/8 Circle: Robotics and Laparoscopic

Needle Points And Cutting Edges










| Needle Type | Point Geometry | Description | Commonly Used For |
|--------------------------------|---|---|--|
| Precision Conventional Cutting |  | Triangular point with cutting edge on the inner curvature and precision-honed for added sharpness | General skin closure, subcutaneous tissue, sometimes for ophthalmic surgery, plastic or reconstructive surgery |
| Reverse Cutting |  | Triangular point with cutting edge on the outer curvature | Skin closure, ligamentous or fibrous tissues |
| Precision Reverse Cutting |  | Triangular point with cutting edge on the outer curvature and precision-honed for added sharpness | Primarily used in plastic surgery for delicate procedures and where a good cosmetic result is desired |
| Taper Point |  | Round shaft, straight or curved, taper point with no cutting edge | Soft tissue closure such as fascia, vascular, gastrointestinal, and most soft tissue below the skin |
| Taper Cutting |  | Four (4) cutting edges with a taper shaft | For penetration through tough tissue |














Precision Conventional Cutting Needles

| Surgical Specialties | Ethicon™ | Covidien™ | |
|----------------------|----------|-----------|---|
| 3/8 Circle | | | |
| PC-1 | PC-1 | PC-13 |  |
| PC-3 | PC-3 | PC-11 |  |
| PC-5 | PC-5 | PC-12 |  |



Precision Reverse Cutting Needles

| Surgical Specialties | Ethicon™ | Covidien™ | |
|----------------------|----------|-----------|---|
| 3/8 Circle | | | |
| P-1 | P-1 | P-10 |  |
| G-1 | G-1 | G-1 |  |
| P-3 | P-3 | P-13 |  |
| G-3 | G-3 | G-3 |  |
| PS-3 | PS-3 | P-11 |  |
| PS-2 | PS-2 | P-12 |  |
| PS-1 | PS-1 | P-14 |  |
| PS | PS | P-15 |  |
| 1/2 Circle | | | |
| PS-4 | PS-4 | P-24 |  |







Reverse Cutting Needles

| Surgical Specialties | Ethicon™ | Covidien™ | |
|----------------------|----------|-----------|---|
| 3/8 Circle | | | |
| C-2 | C-2 | C-2 |  |
| FS-3 | FS-3 | FS-3 |  |
| FS-2 | FS-2 | C-13 |  |
| FS-1 | FS-1 | C-14 |  |
| FS | FS | C-15 |  |
| FSL | FSL | C-16 |  |
| FSLX | FSLX | C-17 |  |
| 1/2 Circle | | | |
| OS-4 | OS-4 | HOS-10 |  |
| CP-2 | CP-2 | GS-10 |  |
| OS-6 | OS-6 | HOS-11 |  |
| CP-1 | CP-1 | GS-11 |  |


















Straight Cutting Needles

| Surgical Specialties | Ethicon™ | Covidien™ | |
|----------------------|----------|-----------|---|
| SC-1 | SC-1 | |  |
| KS | KS | SC-2 |  |

Taper Cutting Needles


| Surgical Specialties | Ethicon™ | Covidien™ | |
|----------------------|----------|-----------|---|
| 3/8 Circle | | | |
| CC | CC | KV-1 |  |
| CC-1 | CC-1 | |  |
| V-4 | V-4 | KV-4 |  |
| 1/2 Circle | | | |
| V-5 | V-5 | KV-5 |  |
| V-7 | V-7 | KV-7 |  |
| V-34 | V-34 | KV-34 |  |

Taper Point Needles



| Surgical Specialties | Ethicon™ | Covidien™ | |
|----------------------|----------|-----------|---|
| 3/8 Circle | | | |
| BV-1 | BV-1 | CV-1 |  |
| BV | BV | CV |  |
| C-1 | C-1 | CV-11 |  |
| BB | BB | CV-15 |  |
| 1/2 Circle | | | |
| RB-2 | RB-2 | CV-22 |  |
| TF | TF | CVF-21 |  |
| RB-1 | RB-1 | CV-23 |  |
| SH-1 | SH-1 | CV-25 |  |
| SH | SH | V-20 |  |
| CT-2 | CT-2 | GS-22 |  |
| MO-6 | MO-6 | HGS-22 |  |
| MH | MH | V-26 |  |
| CT-1 | CT-1 | GS-21 |  |
| CT | CT | GS-24X |  |
| CTX | CTX | GS-25 |  |
| 5/8 Circle | | | |
| UR-6 | UR-6 | GU-46 |  |
| UR-5 | UR-5 | GU-45 |  |

Express Gut™

Precision Conventional Cutting


| | | | Size | | | | | |
|--|--------------|----------|--------|--------|-----|-----|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| PC-1 13mm 3/8 Circle  | Express Gut™ | 18"/45cm | B1916N | B1915N | | | | |

Precision Reverse Cutting






| | | | Size | | | | | |
|---|--------------|----------|--------|--------|-----|-----|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| P-1 11mm 3/8 Circle  | Express Gut™ | 18"/45cm | B1917N | | | | | |
| P-3 13mm 3/8 Circle  | Express Gut™ | 18"/45cm | B1918N | B1919N | | | | |

Plain Gut

Precision Conventional Cutting

| | | | Size | | | | | |
|--|-----------|----------|--------|-----|-----|-----|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| PC-1 13mm 3/8 Circle  | Plain Gut | 18"/45cm | B2711N | | | | | |


Precision Reverse Cutting

| | | | Size | | | | | |
|--|-----------|----------|-------|-------|--------|-------|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| G-1 11mm 3/8 Circle  | Plain Gut | 18"/45cm | B774N | | | | | |
| G-1 11mm 3/8 Circle  | Plain Gut | 18"/45cm | B770N | | | | | |
| P-3 13mm 3/8 Circle  | Plain Gut | 18"/45cm | | B686N | B1644N | | | |
| G-3 13mm 3/8 Circle  | Plain Gut | 18"/45cm | | B772N | | | | |
| PS-2 19mm 3/8 Circle  | Plain Gut | 18"/45cm | | | B1627N | B592N | | |



N=box of 12; D=box of 36

Plain Gut



Reverse Cutting

| | | | Size | | | | | |
|---|-----------|----------|------|-----|-------|-------|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| FS-2 19mm 3/8 Circle  | Plain Gut | 27"/70cm | | | B821N | B822N | | |

Straight Cutting



| | | | Size | | | | | |
|---|-----------|----------|------|-----|--------|-----|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| SC-1 13mm Straight  | Plain Gut | 18"/45cm | | | B1824N | | | |
| SC-1 13mm Straight  | Plain Gut | 18"/45cm | | | B1828N | | | |

Taper Point

| | | | Size | | | | | |
|---|-----------|----------|------|-----|-----|-------|-------|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| CT-1 36mm 1/2 Circle  | Plain Gut | 27"/30cm | | | | B842N | B843N | |
| CT 40mm 1/2 Circle  | Plain Gut | 27"/30cm | | | | B852N | B853N | |








Chromic Gut

| Precision Reverse Cutting  | | | Size | | | | | |
|---|-------------|----------|-------|--------|--------|--------|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| G-1 11mm 3/8 Circle  | Chromic Gut | 12"/30cm | C796N | | | | | |
| G-1 11mm 3/8 Circle  | Chromic Gut | 18"/45cm | C790N | | | | | |
| P-3 13mm 3/8 Circle  | Chromic Gut | 18"/45cm | | C687N | C1654N | | | |
| G-3 13mm 3/8 Circle  | Chromic Gut | 18"/45cm | | C792N | C793N | | | |
| PS-3 16mm 3/8 Circle  | Chromic Gut | 18"/45cm | | C1636N | | | | |
| PS-4 16mm 1/2 Circle  | Chromic Gut | 18"/45cm | | | C1643N | | | |
| PS-2 19mm 3/8 Circle  | Chromic Gut | 18"/45cm | | | C1637N | | | |
| | | 27"/70cm | | | | C1638N | | |

| Reverse Cutting  | | | Size | | | | | |
|---|-------------|----------|------|-----|-------|-------|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| FS-2 19mm 3/8 Circle  | Chromic Gut | 27"/70cm | | | C635N | C636N | | |

N=box of 12; D=box of 36

Chromic Gut





| Taper Point  | | | Size | | | | | |
|--|-------------|----------|-------|-------|-------|-------|-------|-------|
| | Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 1 |
| RB-1  17mm 1/2 Circle | Chromic Gut | 27"/70cm | C202N | C203N | C204N | | | |
| SH-1  22mm 1/2 Circle | Chromic Gut | 27"/70cm | | C181N | C182N | | | |
| SH  26mm 1/2 Circle | Chromic Gut | 27"/70cm | | C121N | C122N | C123N | C124N | |
| CT-1  36mm 1/2 Circle | Chromic Gut | 27"/70cm | | | C810N | C811N | C812N | C813N |
| | | 36"/90cm | | | C922N | C923N | C924N | C925N |
| CT  40mm 1/2 Circle | Chromic Gut | 27"/70cm | | | | | C802N | |
| | | 36"/90cm | | | | C913N | C914N | C915N |
| CTX  48mm 1/2 Circle | Chromic Gut | 36"/90cm | | | | | C904N | C905N |

PolySyn FA™ / Polyglycolic Acid

Precision Reverse Cutting



Size

| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
|--|----------------|----------|-----|-------|-------|-------|-----|---|
| P-1  11mm 3/8 Circle | Undyed Braided | 18"/45cm | | M490N | | | | |
| P-3  13mm 3/8 Circle | Undyed Braided | 18"/45cm | | M493N | M494N | | | |
| PS-2  19mm 3/8 Circle | Undyed Braided | 18"/45cm | | | M496N | M497N | | |
| | | 27"/70cm | | | M426N | | | |
| PS-1  24mm 3/8 Circle | Undyed Braided | 27"/70cm | | | | M442N | | |

Taper Cutting




Size

| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
|---|----------------|----------|-----|-----|-----|-----|-------|---|
| V-34  36mm 1/2 Circle | Undyed Braided | 36"/90cm | | | | | M917N | |




Taper Point






Size

| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
|---|----------------|----------|-----|-----|-----|-------|-------|---|
| CT-1  36mm 1/2 Circle | Undyed Braided | 36"/90cm | | | | M944N | M945N | |


Precision Conventional Cutting 

| | | | Size | | | | | |
|--|---------------------|----------|------|-------|-------|-----|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| PC-1 13mm 3/8 Circle  | Undyed Monofilament | 18"/45cm | | Y834N | | | | |
| PC-3 16mm 3/8 Circle  | Undyed Monofilament | 18"/45cm | | Y844N | Y845N | | | |
| PC-5 19mm 3/8 Circle  | Undyed Monofilament | 18"/45cm | | | Y823N | | | |

Precision Reverse Cutting 

| | | | Size | | | | | |
|--|---------------------|----------|-------|-------|-------|-------|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| P-3 13mm 3/8 Circle  | Undyed Monofilament | 18"/45cm | Y492N | Y493N | Y494N | | | |
| | Violet Monofilament | 18"/45cm | | Y463N | Y464N | | | |
| PS-2 19mm 3/8 Circle  | Undyed Monofilament | 18"/45cm | | Y495N | Y496N | Y497N | | |
| | | 27"/70cm | | | Y426N | Y427N | | |
| | Violet Monofilament | 18"/45cm | | | Y513N | | | |
| PS-1 24mm 3/8 Circle  | Undyed Monofilament | 18"/45cm | | | Y682N | | | |
| | | 27"/70cm | | | Y935N | Y936N | | |





Reverse Cutting 

| | | | Size | | | | |
|--|---------------------|----------|------|-----|-----|-------|---|
| | Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| CP-1 36mm 1/2 Circle  | Undyed Monofilament | 27"/70cm | | | | Y266N | |

N=box of 12; D=box of 36

Monoderm™




Taper Point 

| | | | Size | | | | | |
|---|---------------------|----------|-------|-------|-------|----------------|-------|-------|
| | Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 1 |
| RB-1  17mm 1/2 Circle | Undyed Monofilament | 27"/70cm | Y213D | Y214N | Y215N | | | |
| | Violet Monofilament | 27"/70cm | Y303D | Y304D | Y305D | | | |
| SH  26mm 1/2 Circle | Undyed Monofilament | 27"/70cm | | Y415N | Y416N | Y417N Y417D | | |
| | Violet Monofilament | 27"/70cm | | | Y316N | | | |
| CT-1  36mm 1/2 Circle | Undyed Monofilament | 36"/90cm | | | Y944N | Y945N | Y946N | |
| CTX  48mm 1/2 Circle | Violet Monofilament | 36"/90cm | | | | | Y398N | Y399N |




N=box of 12; D=box of 36

PolySyn™/Polyglycolic Acid

Precision Conventional Cutting

| | | | Size | | | | | |
|--|----------------|----------|------|-------|----------------|-------|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| PC-1 13mm 3/8 Circle  | Undyed Braided | 18"/45cm | | G834N | | | | |
| PC-3 16mm 3/8 Circle  | Undyed Braided | 18"/45cm | | | G845N G835N | | | |
| PC-5 19mm 3/8 Circle  | Undyed Braided | 18"/45cm | | | | G824N | | |






Precision Reverse Cutting

| | | | Size | | | | | |
|--|----------------|----------|-------|-------|--------|-------|-------|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| P-1 11mm 3/8 Circle  | Undyed Braided | 18"/45cm | G489N | | | | | |
| P-3 13mm 3/8 Circle  | Undyed Braided | 18"/45cm | G492N | G493N | G494N | | | |
| | Violet Braided | 18"/45cm | | G463N | G464N | | | |
| PS-3 16mm 3/8 Circle  | Undyed Braided | 18"/45cm | | G500N | | | | |
| PS-4 16mm 1/2 Circle  | Undyed Braided | 18"/45cm | | G506N | G507N | | | |
| PS-2 19mm 3/8 Circle  | Undyed Braided | 18"/45cm | | G495N | G496N | G497N | | |
| | | 27"/70cm | | | G426N | G427N | G428N | |
| | Violet Braided | 27"/70cm | | | G0392N | | | |
| PS-1 24mm 3/8 Circle  | Undyed Braided | 18"/45cm | | | | G683N | | |
| | | 27"/70cm | | | G935N | G936N | | |

N=box of 12; D=box of 36

PolySyn™/Polyglycolic Acid

Reverse Cutting ▾

| | | | Size | | | | | |
|---|----------------|----------|-------|-------|-------|-------|-------|-------|
| | Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 1 |
| FS-2  19mm 3/8 Circle | Undyed Braided | 27"/70cm | G421N | G422N | G423N | | | |
| | Violet Braided | 18"/45cm | | G392N | G393N | | | |
| FS-1  24mm 3/8 Circle | Undyed Braided | 27"/70cm | | G441N | G442N | G443N | | |
| | Violet Braided | 27"/70cm | | G451N | | | | |
| CP-2  26mm 1/2 Circle | Undyed Braided | 27"/70cm | | | | G869N | | |
| CP-1  36mm 1/2 Circle | Undyed Braided | 27"/70cm | | | | G266N | G267N | G268N |
| OS-6  36mm 1/2 Circle | Undyed Braided | 27"/70cm | | | | | | G535N |


N=box of 12; D=box of 36

17surgicalspecialties.com | 877.991.1110

PolySyn™/Polyglycolic Acid







Taper Cutting








| | Material | Length | Size | | | | |
|---|----------------|----------|------|-----|-----|-------|---|
| | | | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| V-34  36mm 1/2 Circle | Undyed Braided | 36"/90cm | | | | G917N | |

Taper Point





| | Material | Length | Size | | | | | |
|---|----------------|----------|-------|-------|-------|-------|-------|---|
| | | | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 1 |
| TF  13mm 1/2 Circle | Undyed Braided | 27"/70cm | G433N | | | | | |
| RB-1  17mm 1/2 Circle | Undyed Braided | 27"/70cm | | G214N | G215N | | | |
| SH  26mm 1/2 Circle | Undyed Braided | 27"/70cm | | G415N | G416N | G417N | | |
| | Violet Braided | 27"/70cm | | | G316N | G317N | | |
| SH-1  22mm 1/2 Circle | Violet Braided | 27"/70cm | | G310N | G311N | | | |
| CT-2  26mm 1/2 Circle | Undyed Braided | 27"/70cm | | | | G269N | G270N | |
| | Violet Braided | 27"/70cm | | | | G333N | G334N | |
| UR-6  26mm 5/8 Circle | Violet Braided | 27"/70cm | | | | G602N | G603N | |

PolySyn™/Polyglycolic Acid



| Taper Point  | | | Size | | | | | |
|--|----------------|----------|------|-----|-------|----------------|-------|-------|
| | Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 1 |
| CT-1  36mm 1/2 Circle | Undyed Braided | 27"/70cm | | | G258N | G259N | G260N | G261N |
| | | 36"/90cm | | | G944N | G945N | G946N | G947N |
| | Violet Braided | 27"/70cm | | | G338N | G339N | G340N | G341N |
| | | 36"/90cm | | | G344N | G345N G354N | G346N | G347N |
| UR-5  36mm 5/8 Circle | Violet Braided | 27"/70cm | | | | | G376N | |
| CT  40mm 1/2 Circle | Undyed Braided | 36"/90cm | | | | G957N | G958N | |
| | Violet Braided | 36"/90cm | | | | | G358N | |
| CTX  48mm 1/2 Circle | Undyed Braided | 36"/90cm | | | | | G978N | G977N |
| | Violet Braided | 36"/90cm | | | | | G370N | G371N |

N=box of 12; D=box of 36




Precision Reverse Cutting 

| | | | Size | | | | | |
|--|---------------------|----------|------|-------|-----------------|-------|-----|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| P-3 13mm 3/8 Circle  | Violet Monofilament | 18"/45cm | | R493N | R494N | | | |
| PS-2 19mm 3/8 Circle  | Violet Monofilament | 18"/45cm | | | R513N/ R496N | R497N | | |







Reverse Cutting 










| | | | Size | | | | | |
|--|---------------------|----------|------|-----|-------|-------|-------|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| FS-2 19mm 3/8 Circle  | Violet Monofilament | 27"/70cm | | | R422N | R423N | | |
| FS-1 24mm 3/8 Circle  | Violet Monofilament | 27"/70cm | | | | | R443N | |






Taper Point 

| | | | Size | | | | |
|--|---------------------|----------|------|-------|-------|-------|-------|
| | Material | Length | 7-0 | 6-0 | 5-0 | 4-0 | 3-0 |
| C-1 13mm 3/8 Circle  | Violet Monofilament | 30"/75cm | | R127N | R126N | | |
| RB-2 13mm 1/2 Circle  | Violet Monofilament | 30"/75cm | | | R148N | | |
| RB-1 17mm 1/2 Circle  | Violet Monofilament | 27"/70cm | | | R303N | R304N | R305N |

Taper Point 

| | | | Size | | | | | |
|--|---------------------|-----------|-------|-----|-------|-------|-------|-------|
| | Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 1 |
| RB-1 17mm 1/2 Circle  | Violet Monofilament | 30"/75cm | R320N | | | | | |
| SH 26mm 1/2 Circle  | Violet Monofilament | 27"/70cm | | | R316N | R317N | | |
| CT-2 26mm 1/2 Circle  | Violet Monofilament | 27"/70cm | | | | R333N | R334N | |
| CT-1 36mm 1/2 Circle  | Violet Monofilament | 27"/70cm | | | | R339N | R340N | R341N |
| | | 36"/90cm | | | | | R346N | R347N |
| CT 40mm 1/2 Circle  | Violet Monofilament | 36"/90cm | | | | | R358N | R359N |
| CTX 48mm 1/2 Circle  | Violet Monofilament | 36"/90cm | | | | | R370N | R371N |
| | | 60"/120cm | | | | | R990N | |


| Precision Reverse Cutting  | | | Size | | | | | |
|---|---------------|----------|--------|-------|-------|--------|--------|-------|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| P-1 11mm 3/8 Circle  | Black Braided | 18"/45cm | D639N | | | | | |
| G-1 11mm 3/8 Circle  | Black Braided | 18"/45cm | D780N | | | | | |
| P-3 13mm 3/8 Circle  | Black Braided | 18"/45cm | D1639N | D640N | D641N | | | |
| G-3 13mm 3/8 Circle  | Black Braided | 18"/45cm | | | D783N | | | |
| PS-1 24mm 3/8 Circle  | Black Braided | 18"/45cm | | | | D1684N | | |
| PS-2 19mm 3/8 Circle  | Black Braided | 18"/45cm | | | D683N | | | |
| PS 26mm 3/8 Circle  | Black Braided | 18"/45cm | | | | | D1588N | |
| PSL 30mm 3/8 Circle  | Black Braided | 18"/45cm | | | | | | D580N |
| | | 30"/75cm | | | | | | D590N |

| Reverse Cutting  | | | Size | | | | | |
|--|---------------|----------|-------|-----|--------|--------|-------|---|
| | Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 1 |
| FS-2 19mm 3/8 Circle  | Black Braided | 18"/45cm | D682N | | D1679N | | | |
| | | 30"/75cm | | | D0768N | | | |
| FS-1 24mm 3/8 Circle  | Black Braided | 18"/45cm | | | D684N | D0238N | | |
| FS 26mm 3/8 Circle  | Black Braided | 18"/45cm | | | | D685N | | |
| | | 30"/75cm | | | D0212N | | | |
| FSL 30mm 3/8 Circle  | Black Braided | 18"/45cm | | | | | D678N | |
| | | 30"/75cm | | | D675N | D679N | D680N | |




N=box of 12; D=box of 36

Silk

Straight Cutting

| | | | Size | | | | | |
|---|---------------|----------|------|-----|-----|-------|----------------|---|
| No Needle | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| KS  60mm Straight | Black Braided | 30"/75cm | | | | D622N | D623N D623D | |

Taper Point

| | | | Size | | | | | |
|---|---------------|----------|-------|-----|--------|-------|-------|-------|
| No Needle | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| C-1  13mm 3/8 Circle | Black Braided | 18"/45cm | D707N | | D7734N | | | |
| SH  26mm 1/2 Circle | Black Braided | 30"/75cm | | | | D832N | D833N | D834N |
| CT-1  36mm 1/2 Circle | Black Braided | 30"/75cm | | | | | | D424N |

Silk Ligatures Tangle Free Dispenser

| | | | Size | | | | | |
|-----------|---------------|-------------|------|-----|-------|-------|-------|-------|
| No Needle | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| | Black Braided | 12x18"/45cm | | | D183D | D184D | D185D | |
| | | 6x30"/75cm | | | | | | D306D |
| | | 12x30"/75cm | | | D303D | D304D | D305D | |

N=box of 12; D=box of 36

Precision Conventional Cutting

| | | | Size | | | | |
|---|----------------------------|----------|--------|--------|--------|--------|-----|
| Needle | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 |
| PC-1 13mm 3/8 Circle  | Black Monofilament | 18"/45cm | A1856N | A1855N | | | |
| | Black Monofilament Pliable | 18"/45cm | | A1955N | | | |
| PC-3 16mm 3/8 Circle  | Black Monofilament | 18"/45cm | | A1865N | A1864N | | |
| | Black Monofilament Pliable | 18"/45cm | | A1965N | | | |
| PC-5 19mm 3/8 Circle  | Black Monofilament | 18"/45cm | | | A1894N | A1893N | |
| | Black Monofilament Pliable | 18"/45cm | | | A1994N | | |




Precision Reverse Cutting

| | | | Size | | | | |
|---|---------------------|----------|--------|-------------------|-----------------|--------|-------|
| Needle | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 |
| P-1 11mm 3/8 Circle  | Black Monofilament | 18"/45cm | A697N | | | | |
| P-3 13mm 3/8 Circle  | Black Monofilament | 18"/45cm | A1698N | A698N AD-698N | A699N | | |
| | Undyed Monofilament | 18"/45cm | | A690N | A691N | | |
| PS-4 16mm 1/2 Circle  | Black Monofilament | 18"/45cm | | | A1662N | | |
| PS-3 16mm 3/8 Circle  | Black Monofilament | 18"/45cm | A1665N | A1668N | | | |
| PS-2 19mm 3/8 Circle  | Black Monofilament | 18"/45cm | A1281N | A1666N AD-1666 | A1667N A662D | A1669N | A593N |
| | Blue Monofilament | 18"/45cm | | | AG667N | | |
| PS-1 24mm 3/8 Circle  | Black Monofilament | 18"/45cm | | | A1629D | A1663N | |








N=box of 12; D=box of 36

Nylon™

Precision Reverse Cutting

| Needle | Material | Length | Size | | | | |
|---|--------------------|----------|------|-----|-----|--------|--------|
| | | | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 |
| PS  26mm 3/8 Circle | Black Monofilament | 18"/45cm | | | | | A585N |
| PSL  30mm 3/8 Circle | Black Monofilament | 30"/75cm | | | | A1691N | |
| PSLX  36mm 3/8 Circle | Black Monofilament | 30"/75cm | | | | A1683N | A1697N |

Reverse Cutting

| | Material | Length | Size | | | | |
|---|--------------------|----------|--------|-------------------|--------|-------------------|--------|
| | | | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 |
| C-2  13mm 3/8 Circle | Black Monofilament | 18"/45cm | A667N | | | | |
| FS-3  16mm 3/8 Circle | Black Monofilament | 18"/45cm | A0181N | A0556N | | | |
| FS-2  19mm 3/8 Circle | Black Monofilament | 18"/45cm | | A661N AC-0568D | A662N | | |
| FS-1  24mm 3/8 Circle | Black Monofilament | 18"/45cm | | | A1629N | A663N | |
| | | 30"/75cm | | | | A669N AC-0586D | |
| FS  26mm 3/8 Circle | Black Monofilament | 18"/45cm | | | | | A664N |
| FSL  30mm 3/8 Circle | Clear Monofilament | 30"/75cm | | | | A1671N | |
| FSLX  36mm 3/8 Circle | Black Monofilament | 30"/75cm | | | | A1673N | A1674N |

N=box of 12; D=box of 36

| Precision Conventional Cutting  | | | | Size | | | | | |
|--|---|-------------------|----------|--------|--------|--------|--------|-----|---|
| | | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| PC-1 13mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | J8617N | J8618N | | | | |
| PC-3 16mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | | J8635N | J8634N | | | |
| PC-5 19mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | | | J8631N | J8632N | | |


| Precision Reverse Cutting  | | | | Size | | | | |
|---|---|-------------------|----------|------|--------|--------|--------|--------|
| | | Material | Length | 7-0 | 6-0 | 5-0 | 4-0 | 3-0 |
| P-1 11mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | | J8697N | | | |
| P-3 13mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | | J8695N | J8698N | J8699N | |
| PS-3 16mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | | J8680N | J8681N | | |
| PS-2 19mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | | | J8686N | J8682N | J8687N |
| PS-1 24mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | | | | | J8663N |

| Reverse Cutting  | | | | Size | | | | |
|---|---|-------------------|----------|------|--------|--------|--------|--------|
| | | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 |
| FS-2 19mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | | J8661N | J8683N | J8665N | |
| FS-1 24mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | | | | J8684N | |
| FS 26mm 3/8 Circle |  | Blue Monofilament | 18"/45cm | | | | | J8685N |





N=box of 12; D=box of 36

Polypropylene


Straight Cutting

| | | | Size | | | | | |
|---|-------------------|----------|------|-----|-----|-----|--------|---|
| | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| KS  60mm Straight | Blue Monofilament | 30"/75cm | | | | | J8623N | |

Taper Cutting

| | | | Size | | | | | |
|--|-------------------|----------|--------|--------|-----|--------|--------|--------|
| | Material | Length | 7-0 | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 |
| CC  9.3mm 3/8 Circle | Blue Monofilament | 24"/60cm | J8704D | J8807D | | | | |
| CC-1  13mm 3/8 Circle | Blue Monofilament | 30"/75cm | | J8707D | | | | |
| V-5  17mm 1/2 Circle | Blue Monofilament | 36"/90cm | | | | J8935D | | |
| V-7  26mm 1/2 Circle | Blue Monofilament | 36"/90cm | | | | J8975D | J8976D | J8977D |

Taper Point







| | | | Size | | | | | |
|--|-------------------|----------|-------------------|--------|-----|-----|-----|-----|
| | Material | Length | 7-0 | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 |
| BV-1  9.3mm 3/8 Circle | Blue Monofilament | 18"/45cm | J8701D | J8806D | | | | |
| | | 24"/60cm | J8702D/ J8704D | J8805D | | | | |
| | | | J8747N | J8807D | | | | |
| | | 30"/75cm | J8703D | J8709D | | | | |
| | | | J8755N | | | | | |

N=box of 12; D=box of 36

27 surgicalspecialties.com | 877.991.1110










Polypropylene

Taper Point 

| | | | Size | | | | | |
|---|-------------------|----------|------|--------|--------|--------|--------|-----|
| | Material | Length | 7-0 | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 |
| C-1  13mm 3/8 Circle | Blue Monofilament | 18"/45cm | | J8718D | J8717D | | | |
| | | 24"/60cm | | J8726D | J8725D | | | |
| | | 30"/75cm | | J8706D | | | | |
| | | 36"/90cm | | | J8720D | | | |
| BV  11mm 3/8 Circle | Blue Monofilament | 24"/60cm | | J8610N | | | | |
| | | 30"/75cm | | J8776N | | | | |
| RB-2  13mm 1/2 Circle | Blue Monofilament | 18"/45cm | | J8714D | J8713D | | | |
| | | 24"/60cm | | | J8554D | | | |
| | | 30"/75cm | | J8711D | J8710D | | | |
| BB  17mm 3/8 Circle | Blue Monofilament | 36"/90cm | | | J8580D | J8581D | | |
| RB-1  17mm 1/2 Circle | Blue Monofilament | 18"/45cm | | | J8756D | | | |
| | | 24"/60cm | | | J8555D | | | |
| | | 36"/90cm | | | J8556D | J8557D | J8558D | |
| SH-1  22mm 1/2 Circle | Blue Monofilament | 30"/75cm | | | | J8526D | J8762D | |

N=box of 12; D=box of 36


Polypropylene

| Taper Point  | | Size | | | | | | | |
|--|-------------------|-----------|--------|--------|--------|--------|--------|--------|--------|
| | | Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 1 |
| SH  26mm 1/2 Circle | Blue Monofilament | 30"/75cm | | | J8832N | J8833N | J8834N | | |
| SH  26mm 1/2 Circle | Blue Monofilament | 36"/90cm | | J8521D | J8522D | J8523D | | | |
| | | 48"/122cm | | | J8534D | J8533D | | | |
| MO-6  26mm 1/2 Circle | Blue Monofilament | 30"/75cm | | | | | | J8418N | |
| CT-2  26mm 1/2 Circle | Blue Monofilament | 30"/75cm | | | | J8411N | J8412N | | |
| MH  36mm 1/2 Circle | Blue Monofilament | 36"/90cm | | | J8842D | J8843D | | | |
| | | 48"/122cm | | | | J8853D | | | |
| CT-1  36mm 1/2 Circle | Blue Monofilament | 30"/75cm | | | | J8423N | J8424N | J8425N | |
| CT  40mm 1/2 Circle | Blue Monofilament | 30"/75cm | | | | | J8434N | J8435N | |
| CTX  48mm 1/2 Circle | Blue Monofilament | 30"/75cm | | | | | | | J8455N |



N=box of 12; D=box of 36

Polyviolene™ XL / Polyester - Coated





Reverse Cutting ▾

| | | Size | | | | | | |
|--|---------------------------|------|-----|-----|-----|---|-------|--|
| Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 2 | |
| OS-4 22mm 1/2 Circle  | Green Braided 30"/75cm | | | | | | E519N | |

Taper Cutting ▾




| | | Size | | | | | | |
|---|---------------------------|--------|-----|-----|-----|-----|--------|---|
| Needle | Material | Length | 6-0 | 5-0 | 4-0 | 3-0 | 2-0 | 0 |
| V-4 17mm 3/8 Circle  | Green Braided 30"/75cm | | | | | | EX927D | |
| V-7 26mm 1/2 Circle  | Green Braided 30"/75cm | | | | | | EX977D | |

Taper Point ⊙

| | | Size | | | | | | |
|--|---------------------------|--------|--------|--------|--------|-------|---|--|
| Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 1 | |
| RB-1 17mm 1/2 Circle  | Green Braided 30"/75cm | | EX551D | | | | | |
| | 36"/90cm | EX556D | EX557D | EX558D | | | | |
| SH-1 22mm 1/2 Circle  | Green Braided 36"/90cm | | | EX762D | EX763D | | | |
| SH 26mm 1/2 Circle  | Green Braided 30"/75cm | | | | EX563D | | | |
| | 36"/90cm | | | | EX523D | E524N | | |
| | White Braided 36"/90cm | | | | EX513D | | | |
| SH 26mm 1/2 Circle  | Green Braided 30"/75cm | | | | E833N | E834N | | |

N=box of 12; D=box of 36

Polyviolene™ XL / Polyester - Coated

| Taper Point  | | | Size | | | | | |
|---|---------------|----------|------|-----|-----|-----|-------|-------|
| | Material | Length | 5-0 | 4-0 | 3-0 | 2-0 | 0 | 1 |
| CT-2  26mm 1/2 Circle | Green Braided | 30"/75cm | | | | | E412N | |
| CT-1  36mm 1/2 Circle | Green Braided | 30"/75cm | | | | | E424N | E425N |

N=box of 12; D=box of 36

Ethicon Cross Reference Guide

| Ethicon | SharpPoint Plus™ | Product Description |
|---------|------------------|--|
| 424H | D424N | Silk 0 Black 30" CT-1 Taperpoint 1/2 Circle 36mm |
| 580H | D580N | Silk 0 Black 18" PSL Reverse Cutting 3/8C 30mm |
| 585H | A585N | Nylon™ 2-0 Black 18" PS Precision Reverse Cutting 3/8C 26mm |
| 590H | D590N | Silk 0 Black 30" PSL Reverse Cutting 3/8C 30mm |
| 593H | A593N | Nylon™ 2-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 622H | D622N | Silk 3-0 Black 30" KS Reverse Cutting Straight 60mm |
| 623H | D623D | Silk 2-0 Black 1x30" KS Straight Cutting 60mm |
| 623H | D623N | Silk 2-0 Black 30" GS60 Straight Straight Cutting 60mm |
| 635H | C635N | Chromic Gut 4-0 27" FS-2 Reverse Cutting 3/8C 19mm |
| 636H | C636N | Chromic Gut 3-0 27" FS-2 Reverse Cutting 3/8C 19mm |
| 639G | D639N | Silk 6-0 Black 18" P-1 Precision Reverse Cutting 3/8C 11mm |
| 639G | D1639N | Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm |
| 640G | D640N | Silk 5-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 641G | D641N | Silk 4-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 661H | A661N | Nylon™ 5-0 Black 18" FS-2 Reverse Cutting 3/8C 19mm |
| 661H | AC0568D | Nylon™ 5-0 Black 18" FS-2 Reverse Cutting 3/8C 19mm |
| 662G | A662N | Nylon™ 4-0 Black 18" FS-2 Reverse Cutting 3/8C 19mm |
| 663G | A663N | Nylon™ 3-0 Black 18" FS-1 Reverse Cutting 3/8C 24mm |
| 664G | A664N | Nylon™ 2-0 Black 18" FS Reverse Cutting 3/8C 26mm |
| 667G | A667N | Nylon 6-0 Black 1x18" C-2 Reverse Cutting 3/8c 13mm |
| 669H | A669N | Nylon™ 3-0 Black 30" FS-1 Reverse Cutting 3/8C 24mm |
| 669H | AC0586D | Nylon™ 3-0 Black 30" FS-1 Reverse Cutting 3/8C 24mm |
| 675H | D675N | Silk 3-0 Black 1x30" FSL Reverse Cutting 3/8c 30mm |
| 678G | D678N | Silk 0 Black 18" FSL Reverse Cutting 3/8C 30mm |
| 679H | D679N | Silk 2-0 Black 30" FSL Reverse Cutting 3/8c 30mm |
| 680H | D680N | Silk 0 Black 30" FSL Reverse Cutting 3/8C 30mm |
| 682G | D682N | Silk, 5-0, Black, 18", FS-2, Reverse Cutting, 3/8 circle, 19mm |
| 684G | D684N | Silk 3-0 Black 18" FS-1 Reverse Cutting 3/8C 24mm |
| 685G | D685N | Silk 2-0 Black 18" FS Reverse Cutting 3/8C 26mm |
| 686G | B686N | Plain Gut 5-0 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 687G | C687N | Chromic Gut 5-0 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 690G | A690N | Nylon™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 691G | A691N | Nylon™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 697H | A697N | Nylon™ 6-0 Black 18" P-1 Precision Reverse Cutting 3/8C 11mm |
| 698G | AD698N | Nylon™ 5-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 698H | A698N | Nylon™ 5-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 699G | A699N | Nylon™ 4-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 707G | D707N | Silk 6-0 Black 18" C-1 Taper Point 13mm 3/8c |
| 770G | B770N | Plain Gut 6-0 18" G-1 Precision Reverse Cutting 3/8C 11mm Double Armed |

Ethicon Cross Reference Guide

| Ethicon | SharpPoint Plus™ | Product Description |
|---------|------------------|--|
| 772G | B772N | Plain Gut 5-0 18" G-3 Precision Reverse Cutting 3/8C 13mm Double Armed |
| 774G | B774N | Plain Gut 6-0 18" G-1 Precision Reverse Cutting 3/8C 11mm |
| 780G | D780N | Silk 6-0 Black 18" G-1 Precision Reverse Cutting 3/8C 11mm Double Armed |
| 783G | D783N | Silk 4-0 Black 18" G-3 Precision Reverse Cutting 3/8C 13mm Double Armed |
| 786G | D1639N | Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm |
| 790G | C790N | Chromic Gut 6-0 18" G-1 Precision Reverse Cutting 3/8C 11mm Double Armed |
| 792G | C792N | Chromic Gut 5-0 18" G-3 Precision Reverse Cutting 3/8C 13mm Double Armed |
| 793G | C793N | Chromic Gut 4-0 18" G-3 Precision Reverse Cutting 3/8C 13mm Double Armed |
| 796G | C796N | Chromic Gut 6-0 18" G-1 Precision Reverse Cutting 3/8C 11mm |
| 802H | C802N | Chromic Gut 0 27" CT Taperpoint 1/2 Circle 40mm |
| 810H | C810N | Chromic Gut 3-0 27" CT-1 Taperpoint 1/2 Circle 36mm |
| 811H | C811N | Chromic Gut 2-0 27" CT-1 Taperpoint 1/2 Circle 36mm |
| 812H | C812N | Chromic Gut 0 27" CT-1 Taperpoint 1/2 Circle 36mm |
| 813H | C813N | Chromic Gut 1 27" CT-1 Taperpoint 1/2 Circle 36mm |
| 842H | B842N | Plain Gut 3-0 27" CT-1 Taperpoint 1/2 Circle 36mm |
| 843H | B843N | Plain Gut 2-0 27" CT-1 Taperpoint 1/2 Circle 36mm |
| 852H | B852N | Plain Gut 3-0 27" CT Taperpoint 1/2 Circle 40mm |
| 853H | B853N | Plain Gut 2-0 27" CT Taperpoint 1/2 Circle 40mm |
| 904H | C904N | Chromic Gut 0 36" CTX Taperpoint 1/2 Circle 48mm |
| 905H | C905N | Chromic Gut 1 36" CTX Taperpoint 1/2 Circle 48mm |
| 913H | C913N | Chromic Gut 2-0 36" CT Taperpoint 1/2 Circle 40mm |
| 914H | C914N | Chromic Gut 0 36" CT Taperpoint 1/2 Circle 40mm |
| 915H | C915N | Chromic Gut 1 36" CT Taperpoint 1/2 Circle 40mm |
| 922H | C922N | Chromic Gut 3-0 36" CT-1 Taperpoint 1/2 Circle 36mm |
| 923H | C923N | Chromic Gut 2-0 36" CT-1 Taperpoint 1/2 Circle 36mm |
| 924H | C924N | Chromic Gut 0 36" CT-1 Taperpoint 1/2 Circle 36mm |
| 925H | C925N | Chromic Gut 1 36" CT-1 Taperpoint 1/2 Circle 36mm |
| 1588H | D1588N | Silk 2-0 Black 18" PS Precision Reverse Cutting 3/8C 26mm |
| 1627H | B1627N | Plain Gut 4-0 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 1629H | A1629N | Nylon™ 4-0 Black 18" FS-1 Reverse Cutting 3/8C 24mm |
| 1636G | C1636N | Chromic Gut 5-0 18" PS-3 Precision Reverse Cutting 3/8C 16mm |
| 1637G | C1637N | Chromic Gut 4-0 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 1638H | C1638N | Chromic Gut 3-0 27" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 1639G | D1639N | Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm |
| 1643G | C1643N | Chromic Gut 4-0 18" PS-4 Precision Reverse Cutting 1/2C 16mm |
| 1644G | B1644N | Plain Gut 4-0 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 1654G | C1654N | Chromic Gut 4-0 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 1662G | A1662N | Nylon™ 4-0 Black 18" PS-4 Precision Reverse Cutting 1/2C 16mm |
| 1663G | A1663N | Nylon™ 3-0 Black 18" PS-1 Precision Reverse Cutting 3/8C 24mm |
| 1665G | A1665N | Nylon™ 6-0 Black 18" PS-3 Precision Reverse Cutting 3/8C 16mm |

Ethicon Cross Reference Guide

| Ethicon | SharpPoint Plus™ | Product Description |
|---------|------------------|--|
| 1666G | A1666N | Nylon™ 5-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 1666G | AD1666 | Nylon™ 5-0 Black 1x18" PS-2 Precision Reverse Cutting 3/8c 19mm |
| 1667G | A1667N | Nylon™ 4-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 1667H | A662D | Nylon™ 4-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 1668G | A1668N | Nylon™ 5-0 Black 18" PS-3 Precision Reverse Cutting 3/8C 16mm |
| 1668G | A0556N | Nylon™ 5-0 Black 1x18" FS-3 Reverse Cutting 3/8c 16mm |
| 1668H | A0556N | Nylon™ 5-0 Black 1x18" FS-3 Reverse Cutting 3/8c 16mm |
| 1669H | A1669N | Nylon™ 3-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 1671H | A1671N | Nylon™ 3-0 Black 30" FSL Reverse Cutting 3/8C 30mm |
| 1673H | A1673N | Nylon™ 3-0 Black 30" FSLX Reverse Cutting 3/8C 36mm |
| 1674H | A1674N | Nylon™ 2-0 Black 30" FSLX Reverse Cutting 3/8C 36mm |
| 1677G | D683N | Silk 4-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 1679H | D1679N | Silk 3-0 Black 1x18" FS-2 Reverse Cutting 3/8c 19mm |
| 1683H | A1683N | Nylon™ 3-0 Black 30" PSLX Reverse Cutting 3/8C 36mm |
| 1690H | D675N | Silk 3-0 Black 30" FSL Precision Reverse Cutting 3/8C 30mm |
| 1691H | A1691N | Nylon™ 3-0 Black 30" PSL Reverse Cutting 3/8C 30mm |
| 1697H | A1697N | Nylon™ 2-0 Black 30" PSLX Reverse Cutting 3/8C 36mm |
| 1698G | A1698N | Nylon™ 6-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 1824H | B1824N | Plain Gut 4-0 18" SC-1 Reverse Cutting Straight 13mm |
| 1828H | B1828N | Plain Gut 4-0 18" SC-1 Reverse Cutting Straight 13mm Double Armed |
| 1855G | A1855N | Nylon™ 5-0 Black 18" PC-1 Precision Conventional Cutting 3/8C 13mm |
| 1856G | A1856N | Nylon™ 6-0 Black 18" PC-1 Precision Conventional Cutting 3/8C 13mm |
| 1864G | A1864N | Nylon™ 4-0 Black 18" PC-3 Precision Conventional Cutting 3/8C 16mm |
| 1865G | A1865N | Nylon™ 5-0 Black 18" PC-3 Precision Conventional Cutting 3/8C 16mm |
| 1893G | A1893N | Nylon™ 3-0 Black 18" PC-5 Precision Conventional Cutting 3/8C 19mm |
| 1894G | A1894N | Nylon™ 4-0 Black 18" PC-5 Precision Conventional Cutting 3/8C 19mm |
| 1915G | B1915N | Express Gut™ 5-0 18" PC-1 Precision Conventional Cutting 3/8C 13mm |
| 1916G | B1916N | Express Gut™ 6-0 18" PC-1 Precision Conventional Cutting 3/8C 13mm |
| 1955G | A1955N | Nylon™ (pliable) 5-0 Black 18" PC-1 Precision Conventional Cutting 3/8C 13mm |
| 1965G | A1965N | Nylon™ (pliable) 5-0 Black 18" PC-3 Precision Conventional Cutting 3/8C 16mm |
| 1994G | A1994N | Nylon™ (Pliable) 4-0 Black 18" PC-5 Precision Conventional Cutting 3/8C 19mm |
| 7734G | D7734N | Silk 4-0 Black 18" C-1 Taper Point 3/8c 13mm |
| 8410H | J8832N | Polypropylene 3-0 Blue 30" SH Taperpoint 1/2 Circle 26mm |
| 8411H | J8833N | Polypropylene 2-0 Blue 30" SH Taperpoint 1/2 Circle 26mm |
| 8423H | J8423N | Polypropylene 2-0 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm |
| 8424H | J8424N | Polypropylene 0 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm |
| 8425H | J8425N | Polypropylene 1 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm |
| 8434H | J8434N | Polypropylene 0 Blue 30" CT Taperpoint 1/2 Circle 40mm |
| 8435H | J8435N | Polypropylene 1 Blue 30" CT Taperpoint 1/2 Circle 40mm |
| 8455H | J8455N | Polypropylene 1 Blue 30" CTX Taperpoint 1/2 Circle 48mm |

Ethicon Cross Reference Guide

| Ethicon | SharpPoint Plus™ | Product Description |
|---------|------------------|--|
| 8521H | J8521D | Polypropylene 4-0 Blue 36" SH Taperpoint 1/2C 26mm Double Armed |
| 8522H | J8522D | Polypropylene 3-0 Blue 36" SH Taperpoint 1/2C 26mm Double Armed |
| 8523H | J8523D | Polypropylene 2-0 Blue 36" SH Taperpoint 1/2C 26mm Double Armed |
| 8526H | J8526D | Polypropylene 4-0 Blue 30" SH-1 Taperpoint 1/2C 22mm Double Armed |
| 8534H | J8534D | Polypropylene 3-0 Blue 48" SH Taperpoint 1/2C 26mm Double Armed |
| 8554H | J8554D | Polypropylene 5-0 Blue 24" RB-2 Taperpoint 1/2C 13mm Double Armed |
| 8555H | J8555D | Polypropylene 5-0 Blue 24" RB-1 Taperpoint 1/2C 17mm Double Armed |
| 8556H | J8556D | Polypropylene 5-0 Blue 36" RB-1 Taperpoint 1/2C 17mm Double Armed |
| 8557H | J8557D | Polypropylene 4-0 Blue 36" RB-1 Taperpoint 1/2C 17mm Double Armed |
| 8558H | J8558D | Polypropylene 3-0 Blue 36" RB-1 Taperpoint 1/2C 17mm Double Armed |
| 8580H | J8580D | Polypropylene 5-0 Blue 36" BB Taperpoint 3/8C 17mm Double Armed |
| 8581H | J8581D | Polypropylene 4-0 Blue 36" BB Taperpoint 3/8C 17mm Double Armed |
| 8617G | J8617N | Polypropylene 6-0 Blue 18" PC-1 Precision Conventional Cutting 3/8C 13mm |
| 8618G | J8618N | Polypropylene 5-0 Blue 18" PC-1 Precision Conventional Cutting 3/8C 13mm |
| 8623H | J8623N | Polypropylene 2-0 Blue 1x30" KS Straight Cutting 60mm |
| 8631G | J8631N | Polypropylene 4-0 Blue 18" PC-5 Precision Conventional Cutting 3/8C 19mm |
| 8632G | J8632N | Polypropylene 3-0 Blue 18" PC-5 Precision Conventional Cutting 3/8C 19mm |
| 8634G | J8634N | Polypropylene 4-0 Blue 18" PC-3 Precision Conventional Cutting 3/8C 16mm |
| 8635G | J8635N | Polypropylene 5-0 Blue 18" PC-3 Precision Conventional Cutting 3/8C 16mm |
| 8661G | J8661N | Polypropylene 5-0 Blue 18" FS-2 Reverse Cutting 3/8C 19mm |
| 8663G | J8663N | Polypropylene 3-0 Blue 18" PS-1 Precision Reverse Cutting 3/8C 24mm |
| 8665G | J8665N | Polypropylene 3-0 Blue 18" FS-2 Reverse Cutting 3/8C 19mm |
| 8680G | J8680N | Polypropylene 6-0 Blue 18" PS-3 Precision Reverse Cutting 3/8C 16mm |
| 8681G | J8681N | Polypropylene 5-0 Blue 18" PS-3 Precision Reverse Cutting 3/8C 16mm |
| 8682G | J8682N | Polypropylene 4-0 Blue 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 8683G | J8683N | Polypropylene 4-0 Blue 18" FS-2 Reverse Cutting 3/8C 19mm |
| 8684G | J8684N | Polypropylene 3-0 Blue 18" FS-1 Reverse Cutting 3/8C 24mm |
| 8685H | J8685N | Polypropylene 2-0 Blue 18" FS Reverse Cutting 3/8C 26mm |
| 8686G | J8686N | Polypropylene 5-0 Blue 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 8687H | J8687N | Polypropylene 3-0 Blue 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| 8695G | J8695N | Polypropylene 6-0 Blue 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 8697G | J8697N | Polypropylene 6-0 Blue 18" P-1 Precision Reverse Cutting 3/8C 11mm |
| 8698G | J8698N | Polypropylene 5-0 Blue 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 8699G | J8699N | Polypropylene 4-0 Blue 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| 8701H | J8701D | Polypropylene 7-0 Blue 18" BV-1 Taperpoint 3/8C 9.3mm Double Armed |
| 8702H | J8702D | Polypropylene 7-0 Blue 24" BV-1 Taperpoint 3/8C 9.3mm Double Armed |
| 8702H | J8747N | Polypropylene 7-0 Blue 24" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| 8703H | J8703D | Polypropylene 7-0 Blue 30" BV-1 Taperpoint 3/8C 9.3mm Double Armed |
| 8704H | J8704D | Polypropylene 7-0 Blue 24" CC Tapercut 3/8C 9.3mm Double Armed |
| 8706H | J8706D | Polypropylene 6-0 Blue 30" C-1 Taperpoint 3/8C 13mm Double Armed |

Ethicon Cross Reference Guide

| Ethicon | SharpPoint Plus™ | Product Description |
|------------|------------------|--|
| 8707H | J8707D | Polypropylene 6-0 Blue 30" CC-1 Tapercut Cardio 3/8C 13mm Double Armed |
| 8709H | J8709D | Polypropylene 6-0 Blue 30" BV-1 Taperpoint 3/8C 9.3mm Double Armed |
| 8710H | J8710D | Polypropylene 5-0 Blue 30" RB-2 Taperpoint 1/2C 13mm Double Armed |
| 8711H | J8711D | Polypropylene 6-0 Blue 30" RB-2 Taperpoint 1/2C 13mm Double Armed |
| 8713H | J8713D | Polypropylene 5-0 Blue 18" RB-2 Taperpoint 1/2C 13mm Double Armed |
| 8714H | J8714D | Polypropylene 6-0 Blue 18" RB-2 Taperpoint 1/2C 13mm Double Armed |
| 8717H | J8717D | Polypropylene 5-0 Blue 18" C-1 Taperpoint 3/8C 13mm Double Armed |
| 8718H | J8718D | Polypropylene 6-0 Blue 18" C-1 Taperpoint 3/8C 13mm Double Armed |
| 8720H | J8720D | Polypropylene 5-0 Blue 36" C-1 Taperpoint 3/8C 13mm Double Armed |
| 8725H | J8725D | Polypropylene 5-0 Blue 24" C-1 Taperpoint 3/8C 13mm Double Armed |
| 8726H | J8726D | Polypropylene 6-0 Blue 24" C-1 Taperpoint 3/8C 13mm Double Armed |
| 8755H | J8755N | Polypropylene 7-0 Blue 30" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| 8756H | J8756D | Polypropylene 5-0 Blue 18" RB-1 Taperpoint 1/2C 17mm Double Armed |
| 8762H | J8762D | Polypropylene 3-0 Blue 30" SH-1 Taperpoint 1/2C 22mm Double Armed |
| 8805H | J8805D | Polypropylene 6-0 Blue 24" BV-1 Taperpoint 3/8C 9.3mm Double Armed |
| 8806H | J8806D | Polypropylene 6-0 Blue 18" BV-1 Taperpoint 3/8C 9.3mm Double Armed |
| 8807H | J8807D | Polypropylene 6-0 Blue 24" CC Tapercut 3/8C 9.3mm Double Armed |
| 8833H | J8411N | Polypropylene 2-0 Blue 30" CT-2 Taperpoint 1/2 Circle 26mm |
| 8834H | J8412N | Polypropylene 0 Blue 30" CT-2 Taperpoint 1/2 Circle 26mm |
| 8834H | J8418N | Polypropylene 0 Blue 30" MO-6 Taperpoint 1/2 Circle 26mm |
| 8834H | J8834N | Polypropylene 0 Blue 30" SH Taperpoint 1/2 Circle 26mm |
| 8842H | J8842D | Polypropylene 3-0 Blue 36" MH Taperpoint 1/2C 36mm Double Armed |
| 8843H | J8843D | Polypropylene 2-0 Blue 36" MH Taperpoint 1/2C 36mm Double Armed |
| 8853H | J8853D | Polypropylene 2-0 Blue 48" MH Taperpoint 1/2C 36mm Double Armed |
| 8935H | J8935D | Polypropylene 4-0 Blue 36" V-5 Tapercut Cardio 1/2C 17mm Double Armed |
| 8975H | J8975D | Polypropylene 4-0 Blue 36" V-7 Tapercut Cardio 1/2C 26mm Double Armed |
| 8976H | J8976D | Polypropylene 3-0 Blue 36" V-7 Tapercut Cardio 1/2C 26mm Double Armed |
| 8977H | J8977D | Polypropylene 2-0 Blue 36" V-7 Tapercut Cardio 1/2C 26mm Double Armed |
| 8886291451 | EX513D | Polyviolene™ XL 2-0 36" SH Taper Point 1/2 Circle 26mm |
| 8886295351 | EX513D | Polyviolene™ XL 2-0 36" SH Taper Point 1/2 Circle 26mm |
| A183H | D183D | Silk 4-0 Black 18" Silk Ligatures No Needle |
| A184H | D184D | Silk 3-0 Black 18" Silk Ligatures No Needle |
| A185H | D185D | Silk 2-0 Black 18" Silk Ligatures No Needle |
| A303H | D303D | Silk 4-0 Black 30" Silk Ligatures No Needle |
| A304H | D304D | Silk 3-0 Black 30" Silk Ligatures No Needle |
| A305H | D305D | Silk 2-0 Black 30" Silk Ligatures No Needle |
| A306H | D306D | Silk 0 Black 30" Silk Ligatures No Needle |
| EP8610H | J8610N | Polypropylene 6-0 Blue 24" BV Taperpoint 3/8 Circle 11mm Double Armed |
| EP8776H | J8776N | Polypropylene 6-0 Blue 30" BV Taperpoint 3/8 Circle 11mm Double Armed |
| G121H | C121N | Chromic Gut 4-0 27" SH Taperpoint 1/2 Circle 26mm |

Ethicon Cross Reference Guide

| Ethicon | Sharpoint Plus™ | Product Description |
|---------|-----------------|--|
| G122H | C122N | Chromic Gut 3-0 27" SH Taperpoint 1/2 Circle 26mm |
| G123H | C123N | Chromic Gut 2-0 27" SH Taperpoint 1/2 Circle 26mm |
| G124H | C124N | Chromic Gut 0 27" SH Taperpoint 1/2 Circle 26mm |
| G181H | C181N | Chromic Gut 4-0 27" SH-1 Taperpoint 1/2 Circle 22mm |
| G182H | C182N | Chromic Gut 3-0 30" HR22 1/2 Circle Taper Point 22mm |
| G667G | AG667N | Nylon™ 4-0 Blue 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| H821H | B821N | Plain Gut 4-0 27" FS-2 Reverse Cutting 3/8C 19mm |
| H822H | B822N | Plain Gut 3-0 27" FS-2 Reverse Cutting 3/8C 19mm |
| J214H | G214N | PolySyn™ 4-0 Undyed 27" RB-1 Taperpoint 1/2 Circle 17mm |
| J215H | G215N | PolySyn™ 3-0 Undyed 27" RB-1 Taperpoint 1/2 Circle 17mm |
| J232H | G416N | PolySyn™ 3-0 Undyed 27" HR26 1/2 Circle Taper Point 26mm |
| J258H | G258N | PolySyn™ 3-0 Undyed 27" CT-1 Taperpoint 1/2 Circle 36mm |
| J259H | G259N | PolySyn™ 2-0 Undyed 27" CT-1 Taperpoint 1/2 Circle 36mm |
| J260H | G260N | PolySyn™ 0 Undyed 27" HR36 1/2 Circle Taper Point 36mm |
| J261H | G261N | PolySyn™ 1 Undyed 27" CT-1 Taperpoint 1/2 Circle 36mm |
| J266H | G266N | PolySyn™ 2-0 Undyed 27" CP-1 Reverse Cutting 1/2C 36mm |
| J267H | G267N | PolySyn™ 0 Undyed 27" CP-1 Reverse Cutting 1/2C 36mm |
| J268H | G268N | PolySyn™ 1 Undyed 27" CP-1 Reverse Cutting 1/2C 36mm |
| J269H | G269N | PolySyn™ 2-0 Undyed 27" HR26 1/2 Circle Taper Point 26mm |
| J269H | G417N | PolySyn™ 2-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm |
| J270H | G270N | PolySyn™ 0 Undyed 27" CT-2 Taperpoint 1/2 Circle 26mm |
| J310H | G310N | Polysyn™ 4-0 Violet 27" SH-1 Taper Point 1/2c 22mm |
| J311H | G311N | Polysyn™ 3-0 Violet 1x27" SH-1 Taper Point 1/2c 22mm |
| J316H | G316N | PolySyn™ 3-0 Violet 27" HR26 1/2 Circle Taper Point 26mm |
| J317H | G317N | PolySyn™ 2-0 Violet 27" HR26 1/2 Circle Taper Point 26mm |
| J333H | G333N | PolySyn™ 2-0 Violet 27" HR26 1/2 Circle Taper Point 26mm |
| J334H | G334N | PolySyn™ 0 Violet 27" CT-2 Taperpoint 1/2 Circle 26mm |
| J338H | G338N | PolySyn™ 3-0 Violet 27" CT-1 Taperpoint 1/2 Circle 36mm |
| J339H | G339N | PolySyn™ 2-0 Violet 27" HR36 1/2 Circle Taper Point 36mm |
| J340H | G340N | PolySyn™ 0 Violet 27" HRX36 1/2 Circle Taper Point 36.6mm |
| J341H | G341N | PolySyn™ 1 Violet 27" HRX36 1/2 Circle Taper Point 36.6mm |
| J344H | G344N | PolySyn™ 3-0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm |
| J345H | G345N | PolySyn™ 2-0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm |
| J345H | G354N | Polysyn™ 2-0 Violet 36" CT-1 Taper Point 1/2c 36mm |
| J346H | G346N | PolySyn™ 0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm |
| J347H | G347N | PolySyn™ 1 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm |
| J358H | G358N | PolySyn™ 0 Violet 36" HRX40 1/2 Circle Taper Point 40mm |
| J370H | G370N | PolySyn™ 0 Violet 36" CTX Taperpoint 1/2 Circle 48mm |
| J371H | G371N | PolySyn™1 Violet 36" CTX Taperpoint 1/2 Circle 48mm |
| J376H | G376N | PolySyn™ 0 Violet 27" UR-5 Taperpoint 5/8C 36mm |

Ethicon Cross Reference Guide

| Ethicon | SharpPoint Plus™ | Product Description |
|---------|------------------|--|
| J392H | G392N | PolySyn™ 4-0 Violet 18" FS-2 Reverse Cutting 3/8C 19mm |
| J393H | G393N | PolySyn™ 3-0 Violet 18" FS-2 Reverse Cutting 3/8C 19mm |
| J415H | G415N | PolySyn™ 4-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm |
| J421H | G421N | Polysyn™ 5-0 Undyed 27" FS-2 Reverse Cutting 3/8c 19mm |
| J422H | G422N | PolySyn™ 4-0 Undyed 27" FS-2 Reverse Cutting 3/8C 19mm |
| J423H | G423N | PolySyn™ 3-0 Undyed 27" FS-2 Reverse Cutting 3/8C 19mm |
| J426H | G426N | PolySyn™ 4-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm |
| J427H | G427N | PolySyn™ 3-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm |
| J428H | G428N | PolySyn™ 2-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm |
| J433H | G433N | Polysyn™ 5-0 Undyed 27" TF Taper Point 1/2c 13mm |
| J441H | G441N | PolySyn™ 4-0 Undyed 27" FS-1 Reverse Cutting 3/8C 24mm |
| J442H | G442N | PolySyn™ 3-0 Undyed 27" FS-1 Reverse Cutting 3/8C 24mm |
| J443H | G443N | PolySyn™ 2-0 Undyed 27" FS-1 Reverse Cutting 3/8C 24mm |
| J451H | G451N | Polysyn™ 4-0 Violet 27" PS-1 Precision Reverse Cutting 3/8c 24mm |
| J463G | G463N | PolySyn™ 5-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| J464G | G464N | PolySyn™ 4-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| J489G | G489N | PolySyn™ 6-0 Undyed 18" P-1 Precision Reverse Cutting 3/8C 11mm |
| J492G | G492N | PolySyn™ 6-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| J493G | G493N | PolySyn™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| J494G | G494N | PolySyn™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| J495G | G495N | PolySyn™ 5-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| J496G | G496N | PolySyn™ 4-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| J497G | G497N | PolySyn™ 3-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| J500G | G500N | PolySyn™ 5-0 Undyed 18" PS-3 Precision Reverse Cutting 3/8C 16mm |
| J506G | G506N | Polysyn, 5-0, Undyed, 18", PS-4, Precision Reverse Cutting, 1/2 Circle, 16mm |
| J507G | G507N | PolySyn™ 4-0 Undyed 18" PS-4 Precision Reverse Cutting 1/2C 16mm |
| J535H | G535N | PolySyn™ 1 Undyed 27" OS-6 Reverse Cutting 1/2C 36mm |
| J602H | G602N | PolySyn™ 2-0 Violet 27" UR-6 Taperpoint 5/8C 26mm |
| J603H | G603N | PolySyn™ 0 Violet 27" UR-6 Taperpoint 5/8C 26mm |
| J662H | G662N | PolySyn™ 4-0 Undyed 27" KS Reverse Cutting Straight 60mm |
| J683G | G683N | PolySyn™ 3-0 Undyed 18" PS-1 Precision Reverse Cutting 3/8C 24mm |
| J824G | G824N | PolySyn™ 3-0 Undyed 18" PC-5 Precision Conventional Cutting 3/8C 19mm |
| J834G | G834N | PolySyn™ 5-0 Undyed 18" PC-1 Precision Conventional Cutting 3/8C 13mm |
| J845G | G835N | Polysyn™ 4-0 Undyed 18" Precision Conventional Cutting 3/8 16mm |
| J845G | G845N | PolySyn™ 4-0 Undyed 18" PC-3 Precision Conventional Cutting 3/8C 16mm |
| J869H | G869N | PolySyn™ 2-0 Undyed 27" CP-2 Reverse Cutting 1/2C 26mm |
| J917H | G917N | PolySyn™ 2-0 Undyed 36" V-34 Tapercut 1/2 Circle 36mm |
| J935H | G935N | PolySyn™ 4-0 Undyed 27" PS-1 Precision Reverse Cutting 3/8C 24mm |
| J936H | G936N | PolySyn™ 3-0 Undyed 27" PS-1 Precision Reverse Cutting 3/8C 24mm |
| J944H | G944N | PolySyn™ 3-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |

Ethicon Cross Reference Guide

| Ethicon | Sharpoint Plus™ | Product Description |
|---------|-----------------|---|
| J945H | G945N | PolySyn™ 2-0 Undyed 36" HR36 1/2 Circle Taper Point 36mm |
| J946H | G946N | PolySyn™ 0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| J947H | G947N | PolySyn™ 1 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| J957H | G957N | PolySyn™ 2-0 Undyed 36" CT Taperpoint 1/2 Circle 40mm |
| J958H | G958N | PolySyn™ 0 Undyed 36" CT Taperpoint 1/2 Circle 40mm |
| J977H | G977N | PolySyn™ 1 Undyed 36" CTX Taperpoint 1/2 Circle 48mm |
| J978H | G978N | PolySyn™ 0 Undyed 36" CTX Taperpoint 1/2 Circle 48mm |
| K832H | D832N | Silk 3-0 Black 30" HR26 1/2 Circle Taper Point 26mm |
| K833H | D833N | Silk 2-0 Black 30" SH Taperpoint 1/2 Circle 26mm |
| K834H | D834N | Silk 0 Black 30" SH Taperpoint 1/2 Circle 26mm |
| U202H | C202N | Chromic Gut 5-0 27" RB-1 Taperpoint 1/2 Circle 17mm |
| U203H | C203N | Chromic Gut 4-0 27" RB-1 Taperpoint 1/2 Circle 17mm |
| U204H | C204N | Chromic Gut 3-0 27" RB-1 Taperpoint 1/2 Circle 17mm |
| VCP944H | M944N | PolySyn FA™ 3-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| VR426 | M426N | PolySyn FA™ 4-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm |
| VR490 | M490N | Polysyn FA™ 5-0 18" P-1 Precision Reverse Cutting 3/8c 11mm |
| VR493 | M493N | PolySyn FA™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| VR494 | M494N | PolySyn FA™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| VR496 | M496N | PolySyn FA™ 4-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| VR497 | M497N | PolySyn FA™ 3-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| VR917 | M917N | Polysyn FA™ 2-0 Undyed 36" V-34 Taper Cutting 1/2c 36mm |
| VR935 | M442N | Polysyn FA™ 3-0 Undyed 27" PS-1 Precision Reverse Cutting 3/8c 24mm |
| VR945 | M945N | PolySyn FA™ 2-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| X424H | E424N | Polyviolene™ 0 Green 30" CT-1 Taperpoint 1/2 Circle 36mm |
| X425H | E425N | Polyviolene™ 1 Green 30" CT-1 Taperpoint 1/2 Circle 36mm |
| X519H | E519N | Polyviolene™ 2 Green 30" OS-4 Reverse Cutting 1/2C 22mm |
| X523H | EX523D | Polyviolene™ 2-0 Green Braided 36" SH Taperpoint 1/2 Circle 26mm |
| X524H | E524N | Polyviolene™ 0 Green 36" SH Taperpoint 1/2 Circle 26mm Double Armed |
| X551H | EX551D | Polyviolene™ 4-0 Green Braided 30" RB-1 Taperpoint 1/2 Circle 17mm |
| X556H | EX556D | Polyviolene™ 5-0 Green Braided 36" RB-1 Taperpoint 1/2 Circle 17mm |
| X557H | EX557D | Polyviolene™ 4-0 Green Braided 36" RB-1 Taperpoint 1/2 Circle 17mm |
| X558H | EX558D | Polyviolene™ 3-0 Green Braided 36" RB-1 Taperpoint 1/2 Circle 17mm |
| X563H | EX563D | Polyviolene™ 2-0 Green Braided 30" SH Taperpoint 1/2 Circle 26mm |
| X762H | EX762D | Polyviolene™ 3-0 Green Braided 36" SH-1 Taperpoint 1/2 Circle 22mm |
| X763H | EX763D | Polyviolene™ 2-0 Green Braided 36" SH-1 Taperpoint 1/2 Circle 22mm |
| X833H | E833N | Polyviolene™ 2-0 Green 30" SH Taperpoint 1/2 Circle 26mm |

Ethicon Cross Reference Guide

| Ethicon | SharpPoint Plus™ | Product Description |
|---------|------------------|--|
| X834H | E412N | Polyviolene™ 0 Green 30" CT-2 Taperpoint 1/2 Circle 26mm |
| X834H | E834N | Polyviolene™ 0 Green 30" SH Taperpoint 1/2 Circle 26mm |
| X927H | EX927D | Polyviolene™ 2-0 Green Braided 30" V-4 Tapercut 3/8 Circle 17mm |
| X977H | EX977D | Polyviolene™ 2-0 Green Braided 30" V-7 Tapercut 1/2 Circle 26mm |
| Y213H | Y213D | Monoderm™ 5-0 Undyed 27" RB-1 Taperpoint 1/2 Circle 17mm |
| Y214H | Y214N | Monoderm™ 4-0 Undyed 27" RB-1 Taperpoint 1/2 Circle 17mm |
| Y215H | Y215N | Monoderm™ 3-0 Undyed 27" RB-1 Taperpoint 1/2 Circle 17mm |
| Y266H | Y266N | Monoderm™ 2-0 Undyed 27" CP-1 Reverse Cutting 1/2C 36mm |
| Y303H | Y303D | Monoderm™ 5-0 Violet 27" RB-1 Taper Point 17mm 1/2C |
| Y304H | Y304D | Monoderm™ 4-0 Violet 27" RB-1 Taper Point 17mm 1/2C |
| Y305H | Y305D | Monoderm™ 3-0 Violet 27" RB-1 Taper Point 17mm 1/2C |
| Y316H | Y316N | Monoderm™ 3-0 Violet 27" SH Taperpoint 1/2 Circle 26mm |
| Y398H | Y398N | Monoderm™ 0 Violet 36" CTX Taperpoint 1/2 Circle 48mm |
| Y399H | Y399N | Monoderm™ 1 Violet 36" CTX Taperpoint 1/2 Circle 48mm |
| Y415H | Y415N | Monoderm™ 4-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm |
| Y416H | Y416N | Monoderm™ 3-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm |
| Y417H | Y417D | Monoderm™ 2-0 Undyed SH Taper Point 1/2c 26mm |
| Y417H | Y417N | Monoderm™ 2-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm |
| Y426H | Y426N | Monoderm™ 4-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm |
| Y427H | Y427N | Monoderm™ 3-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm |
| Y463G | Y463N | Monoderm™ 5-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| Y464G | Y464N | Monoderm™ 4-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| Y492G | Y492N | Monoderm™ 6-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| Y493G | Y493N | Monoderm™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| Y494G | Y494N | Monoderm™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| Y495G | Y495N | Monoderm™ 5-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| Y496G | Y496N | Monoderm™ 4-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| Y497G | Y497N | Monoderm™ 3-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| Y513G | Y513N | Monoderm™ 4-0 Violet 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| Y523H | Y523N | Monoderm™ 3-0 Undyed 27" KS Reverse Cutting Straight 60mm |
| Y682H | Y682N | Monoderm™ 4-0 Undyed 18" PS-1 Precision Reverse Cutting 3/8C 24mm |
| Y823G | Y823N | Monoderm™ 4-0 Undyed 18" PC-5 Precision Conventional Cutting 3/8C 19mm |
| Y834G | Y834N | Monoderm™ 5-0 Undyed 18" PC-1 Precision Conventional Cutting 3/8C 13mm |
| Y844G | Y844N | Monoderm™ 5-0 Undyed 18" PC-3 Precision Conventional Cutting 3/8C 16mm |

Ethicon Cross Reference Guide

| Ethicon | SharpPoint Plus™ | Product Description |
|---------|------------------|--|
| Y845G | Y845N | Monoderm™ 4-0 Undyed 18" PC-3 Precision Conventional Cutting 3/8C 16mm |
| Y935H | Y935N | Monoderm™ 4-0 Undyed 27" PS-1 Precision Reverse Cutting 3/8C 24mm |
| Y936H | Y936N | Monoderm™ 3-0 Undyed 27" PS-1 Precision Reverse Cutting 3/8C 24mm |
| Y944H | Y944N | Monoderm™ 3-0 Undyed 36" CT-1 Taperpoint 1/2 circle 36mm |
| Y945H | Y945N | Monoderm™ 2-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| Y946H | Y946N | Monoderm™ 0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| Z117H | R117N | PDO™ 6-0 Violet 30" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| Z126H | R126N | PDO™ 5-0 Violet 30" C-1 Taperpoint 3/8 Circle 13mm Double Armed |
| Z127H | R127N | PDO™ 6-0 Violet 30" C-1 Taperpoint 3/8 Circle 13mm Double Armed |
| Z135H | R135N | PDO™ 7-0 Violet 30" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| Z148H | R148N | PDO™ 5-0 Violet 30" RB-2 Taperpoint 1/2 Circle 13mm Double Armed |
| Z155H | R155N | PDO™ 7-0 Violet 24" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| Z259H | R259N | PDO™ 2-0 Undyed 27" CT-1 Taperpoint 1/2 Circle 36mm |
| Z303H | R303N | PDO™ 5-0 Violet 27" RB-1 Taperpoint 1/2 Circle 17mm |
| Z304H | R304N | PDO™ 4-0 Violet 27" RB-1 Taperpoint 1/2 Circle 17mm |
| Z305H | R305N | PDO™ 3-0 Violet 27" RB-1 Taperpoint 1/2 Circle 17mm |
| Z316H | R316N | PDO™ 3-0 Violet 27" SH Taperpoint 1/2 Circle 26mm |
| Z320H | R320N | PDO™ 5-0 Violet 30" RB-1 Taperpoint 1/2 Circle 17mm Double Armed |
| Z333H | R333N | PDO™ 2-0 Violet 27" CT-2 Taperpoint 1/2 Circle 26mm |
| Z333H | R317N | PDO™ 2-0 Violet 27" SH Taperpoint 1/2 Circle 26mm |
| Z334H | R334N | PDO™ 0 Violet 27" CT-2 Taperpoint 1/2 Circle 26mm |
| Z339H | R339N | PDO™ 2-0 Violet 27" CT-1 Taperpoint 1/2 Circle 36mm |
| Z340H | R340N | PDO™ 0 Violet 27" CT-1 Taperpoint 1/2 Circle 36mm |
| Z341H | R341N | PDO™ 1 Violet 27" CT-1 Taperpoint 1/2 Circle 36mm |
| Z346H | R346N | PDO™ 0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm |
| Z347H | R347N | PDO™ 1 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm |
| Z358T | R358N | PDO™ 0 Violet 36" CT Taperpoint 1/2 Circle 40mm |
| Z359T | R359N | PDO™ 1 Violet 36" CT Taperpoint 1/2 Circle 40mm |
| Z370T | R370N | PDO™ 0 Violet 36" CTX Taperpoint 1/2 Circle 48mm |
| Z371T | R371N | PDO™ 1 Violet 36" CTX Taperpoint 1/2 Circle 48mm |
| Z416H | R416N | PDO™ 3-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm |
| Z422H | R422N | PDO™ 4-0 Violet 27" FS-2 Reverse Cutting 3/8C 19mm |
| Z423H | R423N | PDO™ 3-0 Violet 27" FS-2 Reverse Cutting 3/8C 19mm |
| Z443H | R443N | PDO™ 2-0 Violet 27" FS-1 Reverse Cutting 3/8C 24mm |
| Z493G | R493N | PDO™ 5-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| Z494G | R494N | PDO™ 4-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm |
| Z496G | R496N | PDO™ 4-0 Violet 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| Z497G | R497N | PDO™ 3-0 Violet 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| Z513G | R513N | PDO™ 4-0 Violet 18" PS-2 Precision Reverse Cutting 3/8C 19mm |
| Z990G | R990N | PDO™ 0 Violet 60" CTX Taperpoint 1/2 Circle 48mm |

Covidien Cross Reference Guide

| Covidien | Sharp Plus™ | Product Description |
|----------|----------------|---|
| CC923 | Y945N | Monoderm™ 2-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| CC924 | Y946N | Monoderm™ 0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| CC975 | Y399N | Monoderm™ 1 Violet 36" CTX Taperpoint 1/2 Circle 48mm |
| CL904 | G370N | PolySyn™ 0 Violet 36" CTX Taperpoint 1/2 Circle 48mm |
| CL905 | G371N | PolySyn™ 1 Violet 36" CTX Taperpoint 1/2 Circle 48mm |
| CL914 | G358N | PolySyn™ 0 Violet 36" CT Taperpoint 1/2 Circle 40mm |
| CL922 | G344N | PolySyn™ 3-0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm |
| CL923 | G345N | PolySyn™ 2-0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm |
| CL924 | G346N | PolySyn™ 0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm |
| CL925 | G347N | PolySyn™ 1 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm |
| CL952 | G944N | PolySyn™ 3-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| CL953 | G945N | PolySyn™ 2-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| CL954 | G946N | PolySyn™ 0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| CL955 | G947N | PolySyn™ 1 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm |
| CL963 | G957N | PolySyn™ 2-0 Undyed 36" CT Taperpoint 1/2 Circle 40mm |
| CL964 | G958N | PolySyn™ 0 Undyed 36" CT Taperpoint 1/2 Circle 40mm |
| CL974 | G978N | PolySyn™ 0 Undyed 36" CTX Taperpoint 1/2 Circle 48mm |
| CL975 | G977N | PolySyn™ 1 Undyed 36" CTX Taperpoint 1/2 Circle 48mm |
| CM904 | Y398N | Monoderm™ 0 Violet 36" CTX Taperpoint 1/2 Circle 48mm |
| CP411 | J8411N | Polypropylene 2-0 Blue 30" CT-2 Taperpoint 1/2 Circle 26mm |
| CP412 | J8412N | Polypropylene 0 Blue 30" CT-2 Taperpoint 1/2 Circle 26mm |
| CP412 | J8418N | Polypropylene 0 Blue 30" MO-6 Taperpoint 1/2 Circle 26mm |
| CP423 | J8423N | Polypropylene 2-0 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm |
| CP424 | J8424N | Polypropylene 0 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm |
| CP425 | J8425N | Polypropylene 1 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm |
| CP434 | J8434N | Polypropylene 0 Blue 30" CT Taperpoint 1/2 Circle 40mm |
| CP435 | J8435N | Polypropylene 1 Blue 30" CT Taperpoint 1/2 Circle 40mm |
| CP455 | J8455N | Polypropylene 1 Blue 30" CTX Taperpoint 1/2 Circle 48mm |
| CS424 | D424N | Silk 0 Black 30" CT-1 Taperpoint 1/2 Circle 36mm |
| GS65M | D834N | Silk 0 Black 30" SH Taperpoint 1/2 Circle 26mm |
| GS66M | D833N | Silk 2-0 Black 30" SH Taperpoint 1/2 Circle 26mm |
| GS67M | D832N | Silk 3-0 Black 30" SH Taperpoint 1/2 Circle 26mm |
| SC5587G | Y463N | Monoderm™ 5-0 Violet 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SC5590G | Y464N | Monoderm™ 4-0 Violet 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SC5616G | Y495N | Monoderm™ 5-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SC5617G | Y496N | Monoderm™ 4-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SC5618G | Y497N | Monoderm™ 3-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |

Covidien Cross Reference Guide

| Covidien | Sharp point Plus™ | Product Description |
|----------|-------------------------|--|
| SC5627G | Y513N | Monoderm™ 4-0 Violet 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SC5688G | Y492N | Monoderm™ 6-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SC5689G | Y493N | Monoderm™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SC5690G | Y494N | Monoderm™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SDN5666G | AG667N | Nylon™ 4-0 Blue 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SL1625G | G489N | PolySyn™ 6-0 Undyed 18" P-1 Precision Reverse Cutting 3/8 Circle 11mm |
| SL1654 | G824N | PolySyn™ 3-0 Undyed 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm |
| SL1657 | G835N | Polysyn™ 4-0 Undyed 18" Precision Conventional Cutting 3/8 16mm |
| SL1657 | G845N | PolySyn™ 4-0 Undyed 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm |
| SL1659 | G843N | Polysyn 5-0 Undyed PC-1 Precision Conventional Cutting 3/8c 13mm |
| SL5626 | G500N | PolySyn™ 5-0 Undyed 18" PS-3 Precision Reverse Cutting 3/8 Circle 16mm |
| SL5627G | G496N | PolySyn™ 4-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SL5628G | G497N | PolySyn™ 3-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SL5632G | G506N | Polysyn 5-0 Undyed 18" PS-4 Precision Reverse Cutting 1/2 Circle 16mm |
| SL5636G | G495N | PolySyn™ 5-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SL5687G | G493N | PolySyn™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SL5688G | G492N | PolySyn™ 6-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SL5691G | G494N | PolySyn™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SL5633G | G507N | PolySyn™ 4-0 Undyed 18" PS-4 Precision Reverse Cutting 1/2 Circle 16mm |
| SL5640G | G683N | PolySyn™ 3-0 Undyed 18" PS-1 Precision Reverse Cutting 3/8 Circle 24mm |
| SN658 | A1629N | Nylon™ 4-0 Black 18" FS-1 Reverse Cutting 3/8 Circle 24mm |
| SN660 | A0181N | Nylon™ 6-0 Black 18" FS-3 Reverse Cutting 3/8c 16mm |
| SN661 | AC-0568D | Nylon™ 5-0 Black 18" FS-2 Reverse Cutting 3/8c 19mm |
| SN661G | A661N | Nylon™ 5-0 Black 18" FS-2 Reverse Cutting 3/8 Circle 19mm |
| SN662G | A662N | Nylon™ 4-0 Black 18" FS-2 Reverse Cutting 3/8 Circle 19mm |
| SN663G | A663N | Nylon™ 3-0 Black 18" FS-1 Reverse Cutting 3/8 Circle 24mm |
| SN669 | A669N | Nylon™ 3-0 Black 30" FS-1 Reverse Cutting 3/8 Circle 24mm |
| SN669 | AC-0586D | Nylon™ 3-0 Black 30" FS-1 Reverse Cutting 3/8c 24mm |
| SN671G | A1671N | Nylon™ 3-0 Black 30" FSL Reverse Cutting 3/8 Circle 30mm |
| SN671G | A1691N | Nylon™ 3-0 Black 30" PSL Reverse Cutting 3/8 Circle 30mm |
| SN764G | A664N | Nylon™ 2-0 Black 18" FS Reverse Cutting 3/8 Circle 26mm |
| SN1693 | A1893N | Nylon™ 3-0 Black 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm |
| SN1964 | A1864N | Nylon™ 4-0 Black 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm |
| SN1994 | A1894N | Nylon™ 4-0 Black 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm |
| SN1994 | A1994N | Nylon™ (Pliable) 4-0 Black 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm |
| SN3697 | A697N | Nylon™ 6-0 Black 18" P-1 Precision Reverse Cutting 3/8 Circle 11mm |
| SN3965 | A1865N | Nylon™ 5-0 Black 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm |
| SN3965 | A1965N | Nylon™ (Pliable) 5-0 Black 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm |

N=box of 12; D=box of 36

Covidien Cross Reference Guide

| Covidien | Sharp point Plus™ | Product Description |
|----------|-------------------------|--|
| SN5662G | A1662N | Nylon™ 4-0 Black 18" PS-4 Precision Reverse Cutting 1/2 Circle 16mm |
| SN5663G | A1663N | Nylon™ 3-0 Black 18" PS-1 Precision Reverse Cutting 3/8 Circle 24mm |
| SN5665G | A1281N | Nylon™ 6-0 Black 18" PS-2 Precision Reverse Cutting 3/8c 19mm |
| SN5666G | A1666N | Nylon™ 5-0 Black 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SN5666G | AD1666 | Nylon™ 5-0 Black 18" PS-2 Precision Reverse Cutting 3/8c 19mm |
| SN5667 | A662D | Nylon™ 4-0 Black 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SN5667G | A1667N | Nylon™ 4-0 Black 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SN5668G | A0556N | Nylon™ 5-0 Black 18" FS-3 Reverse Cutting 3/8c 16mm |
| SN5668G | A1668N | Nylon™ 5-0 Black 18" PS-3 Precision Reverse Cutting 3/8 Circle 16mm |
| SN5668 | A0556N | Nylon™ 5-0 Black 18" FS-3 Reverse Cutting 3/8c 16mm |
| SN5669G | A1669N | Nylon™ 3-0 Black 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SN5670 | A593N | Nylon™ 2-0 Black 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SN5690G | A690N | Nylon™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SN5691 | A691N | Nylon™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SN5696G | A1698N | Nylon™ 6-0 Black 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SN5698G | A698N | Nylon™ 5-0 Black 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SN5698G | AD698N | Nylon™ 5-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm |
| SN5699G | A699N | Nylon™ 4-0 Black 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SP623 | J8623N | Polypropylene 2-0 Blue 30" KS Straight Cutting 60mm |
| SP623BC | J8623N | Polypropylene 2-0 Blue 30" KS Straight Cutting 60mm |
| SP661G | J8661N | Polypropylene 5-0 Blue 18" FS-2 Reverse Cutting 3/8 Circle 19mm |
| SP665G | J8665N | Polypropylene 3-0 Blue 18" FS-2 Reverse Cutting 3/8 Circle 19mm |
| SP683G | J8683N | Polypropylene 4-0 Blue 18" FS-2 Reverse Cutting 3/8 Circle 19mm |
| SP684G | J8684N | Polypropylene 3-0 Blue 18" FS-1 Reverse Cutting 3/8 Circle 24mm |
| SP868G | J8685N | Polypropylene 2-0 Blue 18" FS Reverse Cutting 3/8 Circle 26mm |
| SP1634 | J8634N | Polypropylene 4-0 Blue 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm |
| SP1635 | J8635N | Polypropylene 5-0 Blue 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm |
| SP1643 | J8632N | Polypropylene 3-0 Blue 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm |
| SP1644 | J8631N | Polypropylene 4-0 Blue 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm |
| SP5663G | J8663N | Polypropylene 3-0 Blue 18" PS-1 Precision Reverse Cutting 3/8 Circle 24mm |
| SP5680G | J8680N | Polypropylene 6-0 Blue 18" PS-3 Precision Reverse Cutting 3/8 Circle 16mm |
| SP5681G | J8681N | Polypropylene 5-0 Blue 18" PS-3 Precision Reverse Cutting 3/8 Circle 16mm |
| SP5682G | J8682N | Polypropylene 4-0 Blue 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SP5686G | J8686N | Polypropylene 5-0 Blue 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SP5687G | J8687N | Polypropylene 3-0 Blue 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm |
| SP5695G | J8695N | Polypropylene 6-0 Blue 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |

N=box of 12; D=box of 36

Covidien Cross Reference Guide

| Covidien | Sharpoint Plus™ | Product Description |
|----------|-----------------|--|
| SP5697G | J8697N | Polypropylene 6-0 Blue 18" P-1 Precision Reverse Cutting 3/8 Circle 11mm |
| SP5698G | J8698N | Polypropylene 5-0 Blue 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SP5699 | J8699N | Polypropylene 4-0 Blue 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SP5699G | J8699N | Polypropylene 4-0 Blue 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SS623 | D623D | Silk 2-0 Black 1x30" KS Straight Cutting 60mm |
| SS678 | D580N | Silk 0 Black 18" PSL Reverse Cutting 3/8 Circle 30mm |
| SS678 | D678N | Silk 0 Black 18" FSL Reverse Cutting 3/8 Circle 30mm |
| SS679 | D679N | Silk 2-0 Black 30" FSL Reverse Cutting 3/8c 30mm |
| SS684G | D684N | Silk 3-0 Black 18" FS-1 Reverse Cutting 3/8 Circle 24mm |
| SS685G | D685N | Silk 2-0 Black 18" FS Reverse Cutting 3/8 Circle 26mm |
| SS70M | D590N | Silk 0 Black 30" PSL Reverse Cutting 3/8 Circle 30mm |
| SS70M | D680N | Silk 0 Black 30" FSL Reverse Cutting 3/8 Circle 30mm |
| SS1639G | D639N | Silk 6-0 Black 18" P-1 Precision Reverse Cutting 3/8 Circle 11mm |
| SS1639G | D1639N | Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm |
| SS5639 | D1639N | Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm |
| SS5640 | D640N | Silk 5-0 Black 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SS5641G | D641N | Silk 4-0 Black 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm |
| SS5649G | D1639N | Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm |
| SS5677 | D683N | Silk 4-0 Black 1x18" PS-2 Precision Reverse Cutting 3/8c 19mm |
| SS5679 | D1679N | Silk 3-0 Black 18" FS-2 Reverse Cutting 3/8c 19mm |
| SS5684 | D1684N | Silk 3-0 Black 18" PS-1 Precision Reverse Cutting 3/8 Circle 24mm |
| SS5685G | D0238N | Silk 2-0 Black 18" FS-1 Reverse Cutting 3/8c 24mm |
| SV9915 | M490N | Polysyn™ FA 5-0 18" P-1 Precision Reverse Cutting 3/8c 11mm |
| VP37MX | J8755N | Polypropylene 7-0 Blue 30" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| VP72MX | J8747N | Polypropylene 7-0 Blue 24" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| VP74X | J8610N | Polypropylene 6-0 Blue 24" BV Taperpoint 3/8 Circle 11mm Double Armed |
| VP205X | J8554D | Polypropylene 5-0 Blue 24" RB-2 Taperpoint 1/2 Circle 13mm Double Armed |
| VP422X | J8522D | Polypropylene 3-0 Blue 36" SH Taperpoint 1/2 Circle 26mm Double Armed |
| VP521X | J8521D | Polypropylene 4-0 Blue 36" SH Taperpoint 1/2 Circle 26mm Double Armed |
| VP523 | J8523D | Polypropylene 2-0 Blue 36" SH Taperpoint 1/2 Circle 26mm Double Armed |
| VP542X | J8534D | Polypropylene 3-0 Blue 48" SH Taperpoint 1/2 Circle 26mm Double Armed |
| VP543 | J8533D | Polypropylene 2-0 Blue 48" SH Taperpoint 1/2 Circle 26mm Double Armed |
| VPF543 | J8533D | Polypropylene 2-0 Blue 48" SH Taperpoint 1/2 Circle 26mm Double Armed |
| VP556X | J8556D | Polypropylene 5-0 Blue 36" RB-1 Taperpoint 1/2 Circle 17mm Double Armed |

N=box of 12; D=box of 36

Covidien Cross Reference Guide

| Covidien | Sharp Plus™ | Product Description |
|----------|----------------|--|
| VP557X | J8557D | Polypropylene 4-0 Blue 36" RB-1 Taperpoint 1/2 Circle 17mm Double Armed |
| VP558X | J8558D | Polypropylene 3-0 Blue 36" RB-1 Taperpoint 1/2 Circle 17mm Double Armed |
| VP566X | J8555D | Polypropylene 5-0 24" RB-1 Taper Point 1/2 Circle 17mm |
| VP580X | J8580D | Polypropylene 5-0 Blue 36" BB Taperpoint 3/8 Circle 17mm Double Armed |
| VP581X | J8581D | Polypropylene 4-0 Blue 36" BB Taperpoint 3/8 Circle 17mm Double Armed |
| VP701X | J8701D | Polypropylene 7-0 Blue 18" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| VP702X | J8702D | Polypropylene 7-0 Blue 24" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| VP703X | J8805D | Polypropylene 6-0 Blue 24" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| VP706X | J8706D | Polypropylene 6-0 Blue 30" C-1 Taperpoint 3/8 Circle 13mm Double Armed |
| VP709X | J8709D | Polypropylene 6-0 Blue 30" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| VP710X | J8710D | Polypropylene 5-0 Blue 30" RB-2 Taperpoint 1/2 Circle 13mm Double Armed |
| VP713X | J8703D | Polypropylene 7-0 Blue 30" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| VP714X | J8714D | Polypropylene 6-0 18" RB-2 Taper Point 1/2 Circle 13mm |
| VP715X | J8713D | Polypropylene 5-0 18" RB-2 Taper Point 1/2 Circle 13mm |
| VP718X | J8718D | Polypropylene 6-0 Blue 18" C-1 Taperpoint 3/8 Circle 13mm Double Armed |
| VP719X | J8717D | Polypropylene 5-0 Blue 18" C-1 Taperpoint 3/8 Circle 13mm Double Armed |
| VP720X | J8720D | Polypropylene 5-0 Blue 36" C-1 Taperpoint 3/8 Circle 13mm Double Armed |
| VP725X | J8725D | Polypropylene 5-0 Blue 24" C-1 Taperpoint 3/8 Circle 13mm Double Armed |
| VP726X | J8726D | Polypropylene 6-0 Blue 24" C-1 Taperpoint 3/8 Circle 13mm Double Armed |
| VP733X | J8711D | Polypropylene 6-0 Blue 30" RB-2 Taperpoint 1/2 Circle 13mm Double Armed |
| VP746X | J8806D | Polypropylene 6-0 Blue 18" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed |
| VP760X | J8762D | Polypropylene 3-0 Blue 30" SH-1 Taperpoint 1/2 Circle 22mm Double Armed |
| VP776X | J8776N | Polypropylene 6-0 Blue 30" BV Taperpoint 3/8 Circle 11mm Double Armed |
| VP832X | J8832N | Polypropylene 3-0 Blue 30" SH Taperpoint 1/2 Circle 26mm |
| VP833 | J8833N | Polypropylene 2-0 Blue 30" SH Taperpoint 1/2 Circle 26mm |
| VP834 | J8418N | Polypropylene 0 Blue 30" MO-6 Taperpoint 1/2 Circle 26mm |
| VP834 | J8834N | Polypropylene 0 Blue 30" SH Taperpoint 1/2 Circle 26mm |
| VP842X | J8842D | Polypropylene 3-0 Blue 36" MH Taperpoint 1/2 Circle 36mm Double Armed |
| VP843 | J8843D | Polypropylene 2-0 Blue 36" MH Taperpoint 1/2 Circle 36mm Double Armed |
| VP853 | J8853D | Polypropylene 2-0 48" MH Taper Point 1/2 Circle 36mm |



Order Online www.SurgicalSpecialties.com



Order by Phone 877-991-1110

Sharpoin™ is also available through major healthcare distributors.

**SURGICAL
SPECIALTIES**
CORPORATION™

For More Information, Visit SurgicalSpecialties.com

Surgical Specialties Corporation, 247 Station Drive, Suite NE1, Westwood, MA 02090, U.S.A.
Phone: 610.404.3520 • 877.991.1110 (USA and Canada only) • Fax: 610.404.4010
All trademarks are property of their respective owners. ©2020 Surgical Specialties Corporation.
All Rights Reserved. SSC-028 R8 01/2021